



MRI Guided Electrophysiological Intervention with a Voltage-Based Electro-Anatomic Mapping System

Citation

Tse, Zion, Charles Dumoulin, Ronald Watkins, Israel Byrd, Jeffrey Schweitzer, Raymond Kwong, Gregory F Michaud, and Ehud J Schmidt. 2012. MRI guided electrophysiological intervention with a voltage-based electro-anatomic mapping system. *Journal of Cardiovascular Magnetic Resonance* 14(Suppl 1): P206.

Published Version

doi:10.1186/1532-429X-14-S1-P206

Permanent link

<http://nrs.harvard.edu/urn-3:HUL.InstRepos:9393001>

Terms of Use

This article was downloaded from Harvard University's DASH repository, and is made available under the terms and conditions applicable to Other Posted Material, as set forth at <http://nrs.harvard.edu/urn-3:HUL.InstRepos:dash.current.terms-of-use#LAA>

Share Your Story

The Harvard community has made this article openly available.
Please share how this access benefits you. [Submit a story](#).

[Accessibility](#)



POSTER PRESENTATION

Open Access

MRI guided electrophysiological intervention with a voltage-based electro-anatomic mapping system

Zion Tse^{1*}, Charles Dumoulin², Ronald Watkins³, Israel Byrd⁴, Jeffrey Schweitzer⁴, Raymond Kwong¹, Gregory F Michaud¹, Ehud J Schmidt¹

From 15th Annual SCMR Scientific Sessions
Orlando, FL, USA. 2-5 February 2012

Background

MRI visualizes luminal & vessel-wall anatomy, and identifies edema & scar tissue, contributing to improved electrophysiological (EP) ablative procedures for treatment of Ventricular Tachycardia & Atrial Fibrillation. MRI-guided EP interventions will be performed for the foreseeable future partially in & outside MRI, due to the need for X-ray/Ultrasound-compliant devices. Electromagnetically tracked catheter procedures, today's norm for most EP procedure phases; vascular navigation, Electro-Anatomic-Mapping (EAM, the diagnostic and therapeutic phases), can only be performed outside MRI. Separate MRI tracking is required in MRI, complicating EP procedures which require moving in & out of the bore [1,2]. Continuous catheter tracking using a single system would allow registration-free EAM in & outside MRI. The goal was developing an MR-compatible St. Jude Medical (SJM) EnSite NavX (ESN) voltage-based tracking [3]. ESN applies 5.8/8.0 kHz voltage bursts between 3 pairs of electrodes on the chest, detecting a catheter's position [4], so a challenge for intra-MRI use is MR gradient ramps which interfere with ESN operation. Minimal MR Image Quality (IQ) reduction also needs to be insured, as well as <2oC patient-skin heating due to components in MRI.

Methods

An MR-compatible ESN (Fig.1) minimized electrode heating and IQ reduction inside the scanner room with modified ESN surface-patches, ferrites on coaxial ESN leads (Fig.1(2)) , and RF filters at the penetration panel.

Outside the scanner room, an electronic switching circuit, triggered by a sequence's gradient-ramps and radio-frequency-pulse transmission (GR&RF), disconnected the ESN leads from the ESN receiver when GR&RF was detected, preventing noise from reaching it (Fig.1(3)). 2 SJM MR-compatible deflectable EP catheters were used (Fig.1(4)) [5]. The ESN was tested at 1.5T in heart phantoms and swine models with varying GRE, SSFP & FSE imaging parameters and slice orientations.

Results

Fig.1(5) shows simultaneous ESN & MR tracking of the dual ESN&MR-tracked catheter in the phantom, verifying the lack of interference between the 2 methods. 3 EAM & MRI imaging experiments were performed in swine, also equipped with sensors to measure temperature at the surface ESN electrodes. Fig.1(6) shows EAM during imaging where catheter tracking of 2 catheters (reference catheter in the Coronary Sinus (CS) & mapping catheter in the Left Atrium) was performed simultaneously on the ESN , with <5% positional error of the CS catheter relative to its position outside MRI. IQ Reduction was <5% in both SE & GRE, with tracking updated during 60-80% of sequence duration. TR elongation was required in TR<4ms sequences. Electrode temperature rise was <1oC with 4 Watt/kg SAR sequences. Surface burns did not occur.

Conclusions

An MR-compatible ESN system permits registration-free, minimal-heating, EAM for EP procedure in & outside MRI, with simultaneous imaging possible.

¹Brigham and Women's Hospital, Harvard Medical School, Boston, MA, USA
Full list of author information is available at the end of the article

Fig.1: Overview of MR-compatible ESN System: Inside scanner room - modified surface patches, RF ferrites & low-pass filters. Outside scanner room - GR&RF switching circuit & ESN receiver

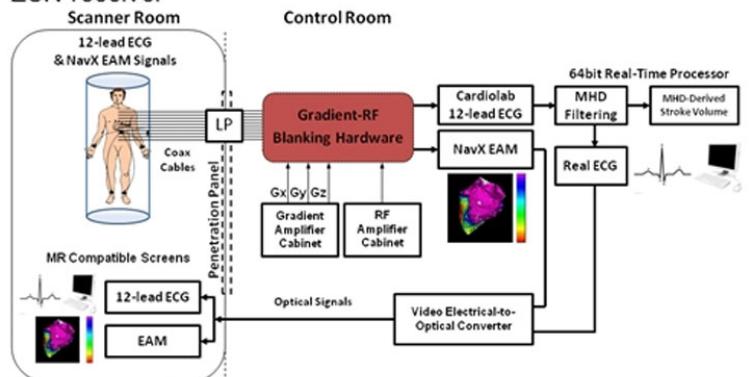


Fig.2: ESN cables with ferrite chucks implemented at RF Max

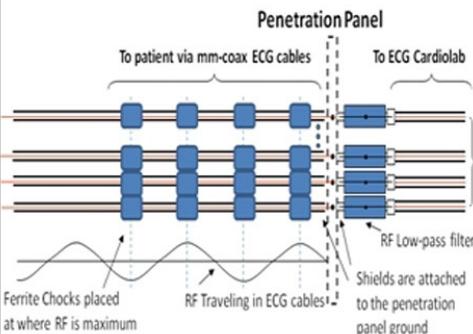


Fig.3: GR-RF blanking circuit for selective EAM signal acquisition (response time<200ns)

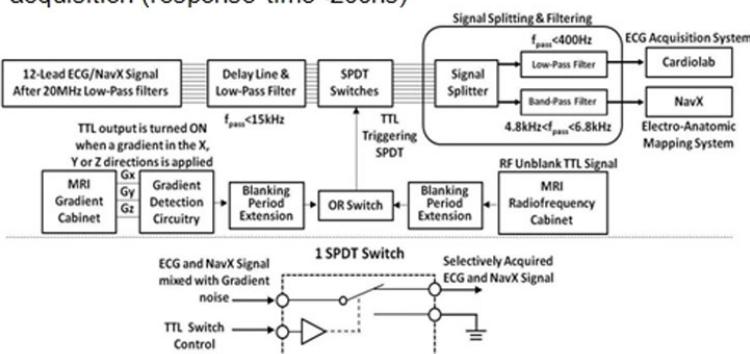


Fig.4: MR-compatible catheters (a) non-irrigated catheters with 5 MR-tracked & 2 NavX electrodes, (b) irrigated with 8 NavX electrodes

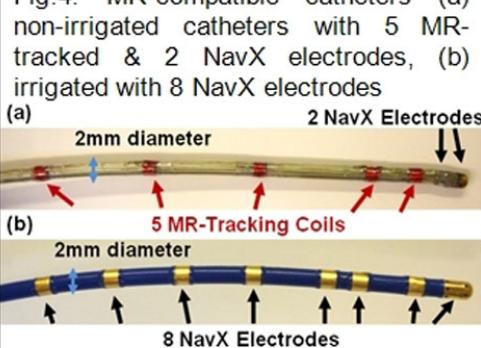


Fig.5: Simultaneous ESN&MR tracking in a heart phantom using an MR- & ESN- tracked catheter

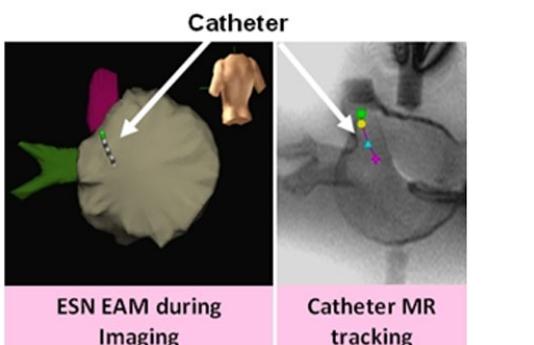


Fig.6: Navx-based EAM during swine MRI imaging. Reference catheter in Coronary Sinus and mapping catheter in Left Atrium

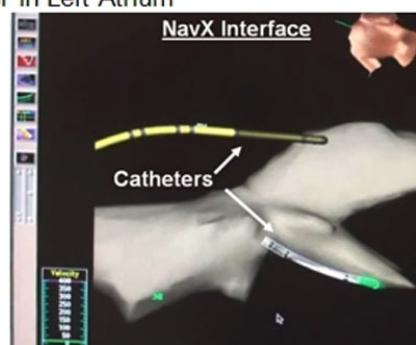


Figure 1

Author details

¹Brigham and Women's Hospital, Harvard Medical School, Boston, MA, USA.
²University of Cincinnati College of Medicine, Cincinnati, OH, USA. ³Stanford University, Stanford, CA, USA. ⁴St Jude Medical, Inc., St Paul, MN, USA.

Published: 1 February 2012

References

1. Dukkipati : Circulation'08..
2. Schmidt : Circulation A&E'09..
3. Nademanee : JACC '04..
4. Early : Eur. Heart J '06..
5. Tse: ISMRM & SCMR '10&'11..

doi:10.1186/1532-429X-14-S1-P206

Cite this article as: Tse et al.: MRI guided electrophysiological intervention with a voltage-based electro-anatomic mapping system. *Journal of Cardiovascular Magnetic Resonance* 2012 14(Suppl 1):P206.