



Improved application of the electrophoretic tissue clearing technology, CLARITY, to intact solid organs including brain, pancreas, liver, kidney, lung, and intestine

Citation

Lee, Hyunsu, Jae-Hyung Park, Incheol Seo, Sun-Hyun Park, and Shin Kim. 2014. "Improved application of the electrophoretic tissue clearing technology, CLARITY, to intact solid organs including brain, pancreas, liver, kidney, lung, and intestine." *BMC Developmental Biology* 14 (1): 781. doi:10.1186/s12861-014-0048-3. <http://dx.doi.org/10.1186/s12861-014-0048-3>.

Published version

<https://doi.org/10.1186/s12861-014-0048-3>

Link

<http://nrs.harvard.edu/urn-3:HUL.InstRepos:13890801>

Terms of use

This article was downloaded from Harvard University's DASH repository, and is made available under the terms and conditions applicable to Other Posted Material (LAA), as set forth at

<https://harvardwiki.atlassian.net/wiki/external/NGY5NDE4ZjgzNTc5NDQzMGIzZWZhMGFIOWI2M2EwYTg>

Accessibility

<https://accessibility.huit.harvard.edu/digital-accessibility-policy>

Share Your Story

The Harvard community has made this article openly available. Please share how this access benefits you. [Submit a story](#)

METHODOLOGY ARTICLE

Open Access

Improved application of the electrophoretic tissue clearing technology, CLARITY, to intact solid organs including brain, pancreas, liver, kidney, lung, and intestine

Hyunsu Lee¹, Jae-Hyung Park², Incheol Seo³, Sun-Hyun Park^{2,5*} and Shin Kim^{4*}

Abstract

Background: Mapping of tissue structure at the cellular, circuit, and organ-wide scale is important for understanding physiological and biological functions. A bio-electrochemical technique known as CLARITY used for three-dimensional anatomical and phenotypical mapping within transparent intact tissues has been recently developed. This method provided a major advance in understanding the structure-function relationships in circuits of the nervous system and organs by using whole-body clearing. Thus, in the present study, we aimed to improve the original CLARITY procedure and developed specific CLARITY protocols for various intact organs.

Results: We determined the optimal conditions for reducing bubble formation, discoloration, and depositing of black particles on the surface of tissue, which allowed production of clearer organ images. We also determined the appropriate replacement cycles of clearing solution for each type of organ, and convincingly demonstrated that 250–280 mA is the ideal range of electrical current for tissue clearing. We then acquired each type of cleared organs including brain, pancreas, liver, lung, kidney, and intestine. Additionally, we determined the images of axon fibers of hippocampal region, the Purkinje layer of cerebellum, and vessels and cellular nuclei of pancreas.

Conclusions: CLARITY is an innovative biochemical technology for the structural and molecular analysis of various types of tissue. We developed improved CLARITY methods for clearing of the brain, pancreas, lung, intestine, liver, and kidney, and identified the appropriate experimental conditions for clearing of each specific tissue type. These optimized methods will be useful for the application of CLARITY to various types of organs.

Keywords: CLARITY, Brain, Nervous system, Electrophoretic tissue clearing, 3D Reconstruction, Purkinje layer

Background

Tissue clearing technologies, such as CLARITY [1] and perfusion-assisted agent release in situ (PARS) [2], which create optically transparent and macromolecule-permeable images, have provided a major advance in the imaging of biological systems. These methods improve tissue permeability by replacing the lipid bilayer of plasma membranes with a nanoporous hydrogel. Unlike mechanical

micro-dissection methods, which aggravate deformation of tissue structure alongside of micro-dissection, CLARITY preserves the intact structure of brain organs, allowing the tracing of neurite projections, and the three-dimensional (3D) and topological reconstruction of traced neurons [1]. PARS can also produce macromolecule permeability and optical transparency in the brain and other organs [2].

Obtaining detailed system-wide informations of organs using a general optical microscope is formidable. For example, single-photon microscopy provides a maximum of 50 μm imaging depth below the organ surface, while even well-optimized two-photon microscopy cannot image deeper than approximately 800 μm [3-5]. However, in various intact organs, CLARITY or PARS can overcome these

* Correspondence: sunhyun_park@hms.harvard.edu; god98005@dsmc.or.kr
Equal contributors

²Department of Physiology, Keimyung University School of Medicine, 1095 Dalgubeoldae-Ro, Dalseo-Gu, Daegu 704-701, South Korea

⁴Department of Immunology, Keimyung University School of Medicine, 1095 Dalgubeoldae-Ro, Dalseo-Gu, Daegu 704-701, South Korea

Full list of author information is available at the end of the article

limitations to permit enhanced viewing of organ structures at greater depths [1,2], providing access to integrated structural and molecular information from the brain and other intact biological systems [1,2]. In particular, CLARITY provides an accelerated rate of un-dissected tissue clearing in the brain through applied electric force, and also avoids tissue damage during clearing [1,2,6].

In the present study, we evaluated whether CLARITY could be effectively applied to other intact organs from the adult mouse (12 weeks-old), including the pancreas, liver, lungs, intestines, and kidneys, as well as the whole brain. Although, this technique has been applied to various cerebral regions of the brain [1], herein, we particularly focused on the Purkinje cell layer of the cerebellum. Furthermore, we examined the structural integrity of cell-vasculature relationships in the pancreatic tail region using single-photon microscopy (0.7 mm). This approach allowed us to obtain a detailed, objective picture of the complexity of the islet vascular system. These data suggest that organ clearing is useful for examining the physiology and pathophysiology of the vascular system as well as other organs.

Methods

Experimental animals

All animal experimental procedures were conducted in accordance with the guidelines of the University Committee on Animal Resources at Keimyung University (Approval No. KM-2014-20R1).

CLARITY applied to the mouse brain

Adult mice (12 weeks old) were anesthetized using a combination of tiletamine-zolazepam-xylazine and perfused transcardially with 30 mL of ice cold 1X phosphate-buffered saline followed by 30 mL of ice cold hydrogel solution with a mixture of 4% PFA, 4% acrylamide, 0.05% bisacrylamide, 0.25% VA044 in PBS. Organs were extracted and incubated in the same solution at 4C for 7 days. The solution temperature was then increased to 37C to initiate polymerization. After 3 h at 37C, hydrogel-embedded organs were placed in an electrophoretic tissue clearing (ETC.) chamber. While sodium borate buffer (200 mM, pH 8.5) containing 4% SDS (the clearing solution) was circulated through the chamber, 250 280 mA was applied across the organs at 42C for 2 4 weeks. The solution circulation velocity was 28 L/min, and the volume of clearing solution was 10 L. The clearing solution replacement cycle is described in Table 1. After clearing, the organs were incubated in PBS at 37C for 4 days to remove SDS. 3-mm-thick horizontal blocks of mouse brain were cleared by electrophoresis for 3 days as described in experimental conditions. The detailed experiment processes were as previously reported [1].

Table 1 Optimized ETC. conditions for various tissues

Organ	Optimized ETC. time (days)	Constant current (mA)	Varying voltage (V)	Replacement interval of clearing solution (days)
Brain	Whole brain: 12-16	280	22 ~ 26	5
	3-mm-thick-block: 3			N/A
Pancreas	8-12	250		7
Kidney	18-22	280		7
Liver	18-22	280		7
Intestine	8-12	250		7
Lung	13-17	250		7

n = 30 for optimization of each organ clearing. ETC., electrophoretic tissue clearing; mA, milliamperes; V, voltage; N/A, not available.

Immunostaining of CLARITY - processed mouse brain and pancreas

To prepare mouse brain for immunostaining, hydrogel-embedded and clarified brains were cut into 3-mm-thick horizontal blocks using a mouse brain matrix (Ted pella). For axon staining, the clarified brain was incubated at 37C for 3 weeks in 0.1% triton X-100, and 1 M sodium borate buffer, pH 8.5, with an anti-tau1 antibody (1:50; Cat# MAB3420, Millipore Corp., Bedford, MA, USA). For pancreatic blood vessels, the clarified pancreas was incubated at 37C for 2 weeks in 0.1% triton X-100 and 1 M sodium borate buffer, pH 8.5, with anti- α -smooth muscle actin (1:50; Cat# ab5694, Abcam, Cambridge, UK). The brain and pancreas were washed at 37C for 1 week in 0.1% triton X-100 and 1 M sodium borate buffer, pH 8.5. The samples were then incubated at 37C for 2 weeks in 0.1% triton X-100 and 1 M sodium borate buffer, pH 8.5, with a goat anti-rabbit-AlexaFluor488 secondary antibody (1:100; Cat# A11029, Invitrogen). For nuclear staining, the clarified pancreas was incubated at 37C for 30 min in DAPI (0.1 μ g/mL, Cat# D9542, Sigma-Aldrich, St. Louis, MO, USA).

Imaging of CLARITY - processed mouse brain and pancreas

For imaging the structures of brain and pancreas, the clarified brain and pancreas were incubated in FocusClear, a water-based immersion medium, for 4 days. The brain and pancreas were then enclosed between two coverglass-bottom petri dishes. The brain and pancreas were imaged (Z-stack volume, 110 650 μ m) using confocal laser scanning microscopy (LSM 5 EXCITER; Carl Zeiss, Jena, Germany) at excitation wavelengths of 488 nm with a C-Apochromat 40 objective (1.2 numerical aperture). Three-dimensional reconstruction was generated using LSM5 EXCITER software (Carl Zeiss). To extract the pancreatic microvasculature images, we carefully traced the anti- α -smooth muscle actin profiles using the ImageJ plug-in Single Neurite Tracer [7].

Results

Electrical current conditions related to enhancing CLARITY

To clear the hydrogel-embedded tissues after completing the polymerization, an electric field was applied to the tissues in an ETC chamber which was built as described previously [1]. Although the previous study recommended using high voltage and temperature to clear tissue, the high constant voltage caused bubble formation, discoloration and/or the deposition of black particles in tissue samples (Additional file 1). Specifically, the

constant voltage created a high current thus damaging the tissue. The previous study also recommended low voltage (V) (10 - 40 V) for ETC. [1], but this condition was shown to be ineffective in applying to other types of tissue as well as the whole brain. Therefore, we determined the optimal range of electrical current condition for clearing tissues.

We found that for brain samples, 280 mA was the optimal current without causing any damage. When the constant current of 280 mA was maintained across the

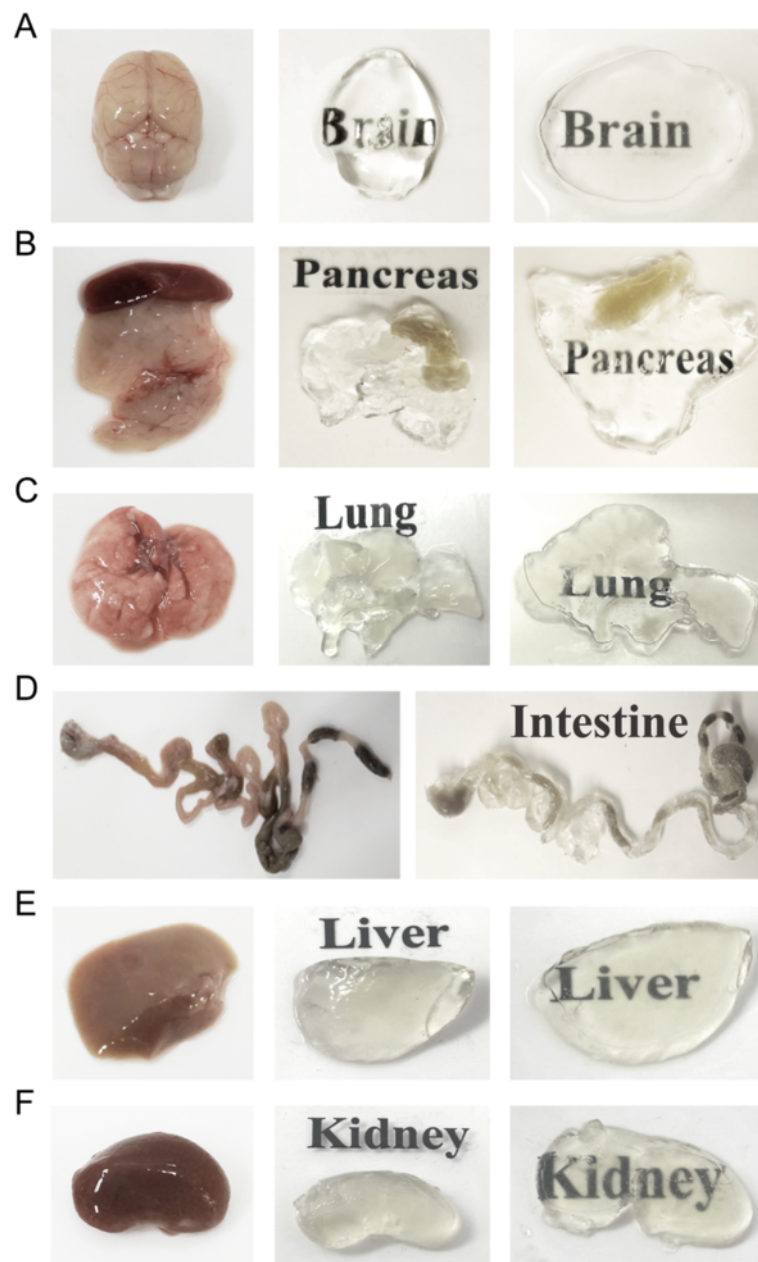


Figure 1 CLARITY of intact adult mouse tissues. In adult mouse tissues (12 weeks old), imaging was performed before and/or after CLARITY. (A) Brain. (B) Pancreas. (C) Lung. (D) Intestine. (E) Liver. (F) Kidney.

sample and the clearing solution was replaced three times over 2 weeks, the sample was completely cleared (Figure 1A). Additionally, it was determined that at a constant electrical current condition, the measured voltage between the platinum wires of chamber should reach at least 20–30 V for successful tissue clearing.

We then investigated the optimal electrical current conditions for the clearing of mouse pancreas and lung tissues. When we identified the appropriate clearing conditions of the tissues, 250 mA is the appropriate current because these tissues are softer than brain tissue. The pancreas sample was cleared by applying 250 mA and replacing the clearing solution two times over 2 weeks (Figure 1B). Lung tissue was cleared by applying 250 mA and replacing the clearing solution three times over 3 weeks (Figure 1C). We tested additional solid organs for clearing under these conditions. The mouse intestine sample was cleared by applying 250 mA and replacing the clearing solution three times over 3 weeks (Figure 1D). Finally, mouse liver and kidney samples were cleared by applying 280 mA and replacing the clearing solution four times over 4 weeks (Figure 1E, F).

Imaging of structures in the adult mouse brain and intact pancreas

We performed immunostaining in adult mouse brain using a microtubule-associated protein tau-antibody. Using confocal images, tau-stained axon fiber was confirmed in the hippocampal region (Additional file 2) and the Purkinje layer of cerebellum (Figure 2). As shown in Additional file 2, the hippocampal region was reconstructed in 3D stacks of images (Z-stack volume, 650 μm). Moreover, the Purkinje layer of cerebellum was reconstructed in 3D stacks of images (Z-stack volume, 210 μm) (Additional file 3 and Additional file 4). Furthermore, we performed imaging in pancreas sample with an α -smooth muscle actin-antibody which is the marker for blood vessels including capillary vessels [8]. We could identify stained vessels in randomly selected regions of distal pancreas in optically transparent pancreas (Figure 3A, Additional files 5, 6, 7, 8 and 9). The vessel and pancreatic region was reconstructed in 3D stacks of images (Z-stack volume, 650 μm). To reconstruct the 3D morphology of the pancreatic capillaries, we traced the profile of the α -smooth muscle actin and extracted the image of the pancreatic capillary using ImageJ (Figure 3B, Additional file 10). To identify the structural integrity of cells-vasculature relationships in distal region of the pancreas, we fluorescently stained the cleared pancreas samples with an α -smooth muscle actin-antibody and DAPI, which binds strongly to A-T rich regions in DNA. We could identify stained vessels and pancreatic cellular nuclei in randomly selected regions of the distal pancreas in optically transparent pancreas (Z-stack volume, 110 μm) (Figure 4, Additional file 11).

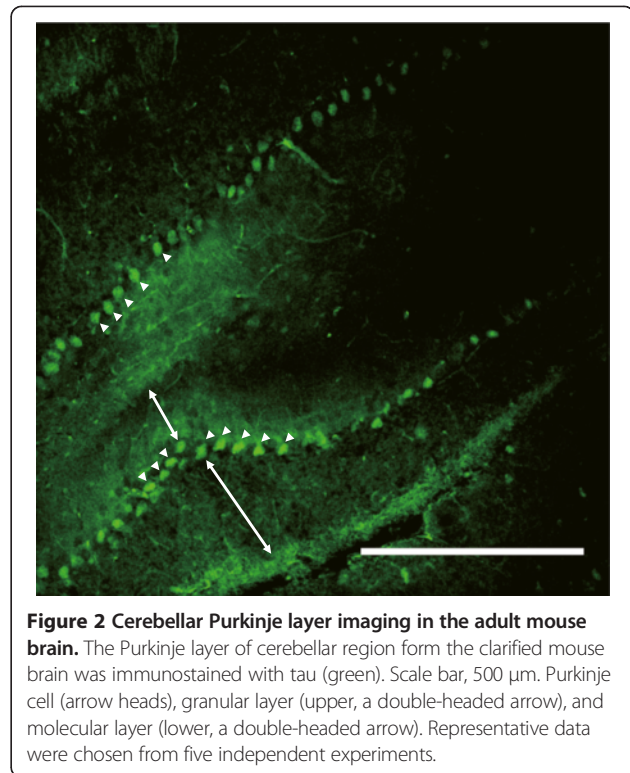


Figure 2 Cerebellar Purkinje layer imaging in the adult mouse brain. The Purkinje layer of cerebellar region from the clarified mouse brain was immunostained with tau (green). Scale bar, 500 μm . Purkinje cell (arrow heads), granular layer (upper, a double-headed arrow), and molecular layer (lower, a double-headed arrow). Representative data were chosen from five independent experiments.

Discussion

Using our optimized CLARITY protocol, we were able to successfully image various mouse whole organs, including the brain and the pancreas. Imaging of blood vessels in the pancreas allowed the general vessel pattern to be determined under normal conditions. This clearing technique may be useful in future studies related to pancreatic disease, of which there have been limited images of blood vessels or other structures. We were also able to identify tau-stained neurons in the Purkinje layer of the cerebellum. Further, we determined the optimal electrical current conditions for clearing of specific organs without bubbling, discoloration, or accumulation of black particles in tissue samples during ETC.

A previous study reported that CLARITY is an effective imaging technique for mapping the nervous system in the brain [1]. However, as bubbling, discoloration, and the appearance of black particles were observed in some samples, further investigation to determine optimal clearing conditions were recommended [1]. In the present study, we found that the range of voltage (10–60 V) was too broad to be practical. Thus, we determined the appropriate experimental conditions for a variety of tissues to produce more clear images and minimize tissue damage, and found an optimal electrical current of 250–280 mA (Table 1). Yang et al. reported that PARS method, a method that provides tissue clearance of intact whole-

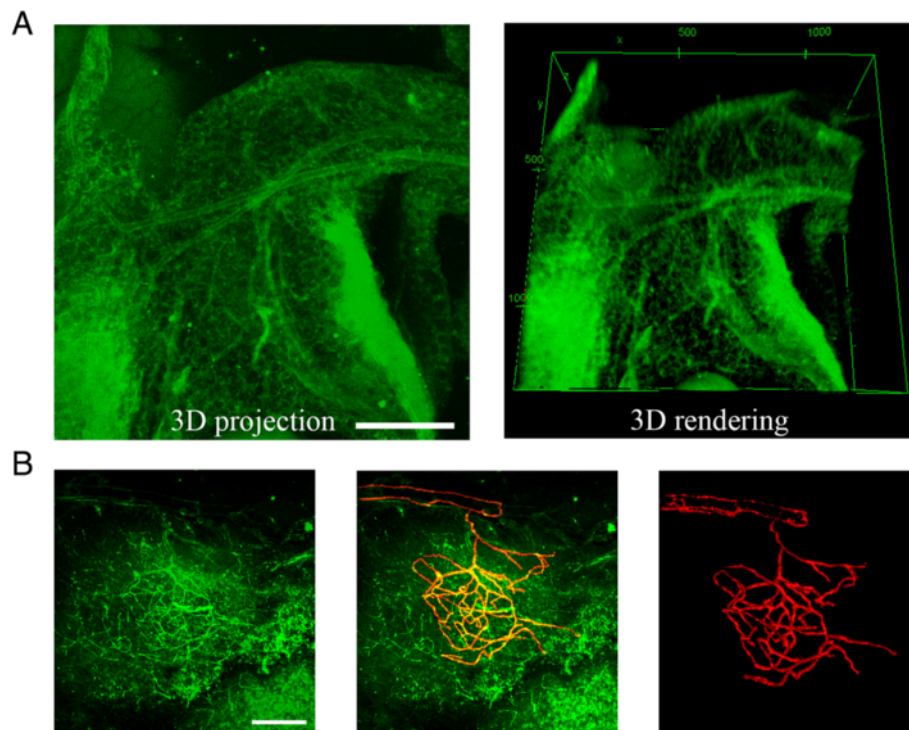


Figure 3 Pancreatic vessel imaging in the intact adult mouse pancreas. In adult mouse tissues (12 weeks old), imaging was performed after CLARITY. **(A)** Three-dimensional (3D) projection (left panel) and rendering (right panel) of clarified mouse pancreas without capillary immunostained for α -smooth muscle actin (green). Scale bar, 300 μ m (Additional files 3 and 4). **(B)** Three-dimensional projection (left panel) clarified mouse pancreas with capillary immunostained for α -smooth muscle actin (green). Scale bar, 200 μ m. Merged image with manually traced (green, middle panel) and 3D rendering extracted capillary (right panel). Representative data were chosen from five independent experiments.

organisms, could provide access to integrated structural and molecular information from the brain and other intact biological systems [2]. Although, organs with a higher cell density (e.g., kidney and liver) required up to 22 days for tissue clearance, using our appropriate experimental conditions, we demonstrated that CLARITY with ETC. could be effectively applied to organs such as the pancreas, liver, lung, kidney and intestine without dissecting the organs. Although Yang et al. introduced a protocol for passive tissue clearing passive clarity technique (PACT), the method

required 3 weeks for clearing of un-dissected, whole adult mouse brain without electrophoretic instrumentation (passive CLARITY) [2]. In our study, although longer times were required to clear tissues than those previously reported [1], our methods provide an important advance for the CLARITY technique, allowing application to a variety of other types of organ.

A major finding of our study was that CLARITY could be successfully applied to intact mouse organs including the brain, pancreas, liver, lung, kidney, and intestine

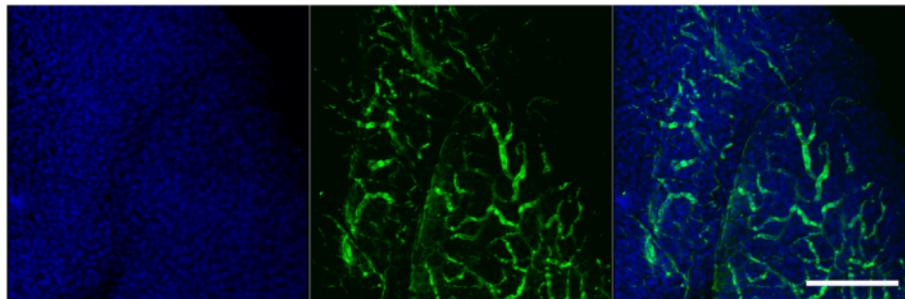


Figure 4 Pancreatic vessel imaging in the intact adult mouse pancreas. Three-dimensional projection (left panel and middle panel) clarified mouse pancreas with nuclear stained for DAPI (blue) and vascular immunostained for α -smooth muscle actin (green). Merged image with 3D rendering nuclei (blue, left panel) and 3D rendering vessels (green, middle panel). Scale bar, 200 μ m. Representative data were chosen from five independent experiments.

(Figure 1). For example, we were able to visualize intact structures and reconstruct the pattern of blood vessels with 3D reconstruction in the optically transparent pancreas. The vascular pattern of the intact pancreas is poorly understood as it is challenging to trace whole vessels using general imaging with 3D reconstruction. Thus, our modified technique may be particularly useful for studying the vascular system in pancreatic function and homeostasis [4,9]. We also demonstrated the 3D reconstruction of the Purkinje layer and adjacent layers in the cerebellar region, suggesting that this modified technique can be applied to the understanding of human disease related to Purkinje cell.

In the present study, we have optimized the CLARITY method for use in tissue clearing of the pancreas, liver, lung, kidney and intestine. In particular, we confirmed that the clarified pancreas allows imaging of the partial vascular structure using single-photon microscopy. Further, we provide the first evidence of 3D imaging of the Purkinje layer of cerebellar region and the cell-vasculature relationships. Overall, our data suggest that CLARITY may be useful for examining a variety of biological tissues.

Conclusions

The technique known as CLARITY used for 3D anatomical and phenotypical mapping within transparent intact tissues has been recently developed. This enhanced technology provides an innovative method for the acquisition of integrated structural and molecular information from the various intact organs across species. The results of this study show that the applicability of CLARITY with ETC. to other types of organ including the pancreas, lung, intestine, liver, and kidney. During establishing the suitable protocols for each intact organ, we encountered difficult situations such as bubble formation, discoloration, and the depositing of block particles on the surface of tissue. However, we overcame the experimental troubles and elaborated the optimal conditions, thus producing more clear images. Importantly, we determined that the supply of 250–280 mA is the ideal range of electrical current for tissue clearing. Moreover, this study demonstrated for the first time that the Purkinje layer of cerebellar region and the cells-vasculature relationships of pancreatic region were reconstructed in 3D stacks of image (Z-stack volume, 650 μm). Finally, having obtained the desired image using CLARITY, we concluded that this method can be applied to various organs besides mouse brain and it would be helpful for other researchers who want to conduct their own CLARITY research.

Additional files

Additional file 1: Figure that Problems of high voltage condition (30 voltage, 42C) during CLARITY. (A) Bubble formation and discoloration in brain. (B) Deposit of black particle in brain.

Additional file 2: Figure that three-dimensional (3D) projection (left panel) and rendering (right panel) of clarified mouse brain immunostained for tau (green). Scale bar, 300 μm .

Additional file 3: Movie that 3D immunohistological visualization of the Purkinje layer of cerebellar region in Figure 2.

Additional file 4: Movie that 3D immunohistological visualization of the Purkinje layer of cerebellar region in Figure 2.

Additional file 5: Movie that 3D immunohistological visualization of alpha-smooth muscle actin in mouse pancreas.

Additional file 6: Movie that 3D immunohistological visualization of alpha-smooth muscle actin in mouse pancreas.

Additional file 7: Movie that 3D immunohistological visualization of alpha-smooth muscle actin in mouse pancreas.

Additional file 8: Movie that 3D immunohistological visualization of alpha-smooth muscle actin in mouse pancreas.

Additional file 9: Movie that 3D immunohistological visualization of alpha-smooth muscle actin in mouse pancreas.

Additional file 10: Movie that 3D visualization of extracted capillary in Figure 3B.

Additional file 11: Movie that 3D visualization of vessels and extracted capillary in Figure 4.

Abbreviations

mA: milliampere; ETC.: Electrophoretic tissue clearing; PARS: Perfusion-assisted agent release in situ; PACT: Passive clarity technique.

Competing interests

The authors declare that they have no competing interest.

Authors contributions

HL, JHP contributed to the conception and design of the study, acquisition of data, interpretation of results and the writing of the manuscript. ICS contributed to the acquisition of data. SHP and SK contributed to the conception of the study, the interpretation of results and writing of the manuscript. All authors gave final approval of the version to be published.

Acknowledgements

We thank all our colleagues in Keimyung University School of Medicine for helpful advice and discussions. This work was supported by the National Research Foundation of Korea (NRF) Grant funded by the Korean Government (MSIP) (No. 2014R1A5A2010008).

Duality of interest

The authors declare that there is no duality of interest associated with this manuscript.

Author details

¹Department of Anatomy, Keimyung University School of Medicine, 1095 Dalgubeoldae-Ro, Dalseo-Gu, Daegu 704-701, South Korea. ²Department of Physiology, Keimyung University School of Medicine, 1095 Dalgubeoldae-Ro, Dalseo-Gu, Daegu 704-701, South Korea. ³Department of Microbiology, Keimyung University School of Medicine, 1095 Dalgubeoldae-Ro, Dalseo-Gu, Daegu 704-701, South Korea. ⁴Department of Immunology, Keimyung University School of Medicine, 1095 Dalgubeoldae-Ro, Dalseo-Gu, Daegu 704-701, South Korea. ⁵Department of Microbiology and Immunobiology, Harvard Medical School, Division of Neuroscience New England Primate Research Center, 1 Pine Hill Drive, Southborough, MA 01772, USA.

Received: 1 August 2014 Accepted: 11 December 2014

Published online: 21 December 2014

References

1. Chung K, Wallace J, Kim SY, Kalyanasundaram S, Andalman AS, Davidson TJ, Mirzabekov JJ, Zalocusky KA, Mattis J, Denisin AK, Pak S, Bernstein H, Ramakrishnan C, Grosenick L, Gradinaru V, Deisseroth K: **Structural and molecular interrogation of intact biological systems.** *Nature* 2013, 497:332–337.

2. Yang B, Treweek JB, Kulkarni RP, Deverman BE, Chen CK, Lubeck E, Shah S, Cai L, Gradinaru V: **Single-cell phenotyping within transparent intact tissue through whole-body clearing.** *Cell* 2014, **158**:945-958.
3. Helmchen F, Denk W: **Deep tissue two-photon microscopy.** *Nat Methods* 2005, **2**:932-940.
4. Rodriguez-Diaz R, Abdulreda MH, Formoso AL, Gans I, Ricordi C, Berggren PO, Caicedo A: **Innervation patterns of autonomic axons in the human endocrine pancreas.** *Cell Metab* 2011, **14**:45-54.
5. Speier S, Nyqvist D, Kohler M, Caicedo A, Leibiger IB, Berggren PO: **Noninvasive high-resolution in vivo imaging of cell biology in the anterior chamber of the mouse eye.** *Nat Protoc* 2008, **3**:1278-1286.
6. Tomer R, Ye L, Hsueh B, Deisseroth K: **Advanced CLARITY for rapid and high-resolution imaging of intact tissues.** *Nat Protoc* 2014, **9**:1682-1697.
7. Abramoff MD, Magalhães PJ, Ram SJ: **Image processing with ImageJ.** *Biophoton Int* 2004, **11**:36-43.
8. Skalli O, Pelte MF, Pecllet MC, Gabbiani G, Gugliotta P, Bussolati G, Ravazzola M, Orci L: **Alpha-smooth muscle actin, a differentiation marker of smooth muscle cells, is present in microfilamentous bundles of pericytes.** *J Histochem Cytochem* 1989, **37**:315-321.
9. Speier S, Nyqvist D, Cabrera O, Yu J, Molano RD, Pileggi A, Moede T, Kohler M, Wilbertz J, Leibiger B, Ricordi C, Leibiger IB, Caicedo A, Berggren PO: **Noninvasive in vivo imaging of pancreatic islet cell biology.** *Nat Med* 2008, **14**:574-578.

Submit your next manuscript to BioMed Central and take full advantage of:

- Convenient online submission
- Thorough peer review
- No space constraints or color figure charges
- Immediate publication on acceptance
- Inclusion in PubMed, CAS, Scopus and Google Scholar
- Research which is freely available for redistribution

Submit your manuscript at
www.biomedcentral.com/submit

