



# Imaging surrogate of the ventricular arrhythmia substrate in a swine model of ventricular tachycardia: high resolution LGE vs. high resolution electro-anatomical mapping

## Citation

Roujol, Sébastien, Tamer A Basha, Cory M Tschabrunn, Kraig V Kissinger, Mark E Josephson, Warren J Manning, Elad Anter, and Reza Nezafat. 2015. "Imaging surrogate of the ventricular arrhythmia substrate in a swine model of ventricular tachycardia: high resolution LGE vs. high resolution electro-anatomical mapping." *Journal of Cardiovascular Magnetic Resonance* 17 (1): Q133. doi:10.1186/1532-429X-17-S1-Q133. <http://dx.doi.org/10.1186/1532-429X-17-S1-Q133>.

## Published version

<https://doi.org/10.1186/1532-429X-17-S1-Q133>

## Link

<http://nrs.harvard.edu/urn-3:HUL.InstRepos:14065522>

## Terms of use

This article was downloaded from Harvard University's DASH repository, and is made available under the terms and conditions applicable to Other Posted Material (LAA), as set forth at

<https://harvardwiki.atlassian.net/wiki/external/NGY5NDE4ZjgzNTc5NDQzMGIzZWZhMGFIOWI2M2EwYTg>

## Accessibility

<https://accessibility.huit.harvard.edu/digital-accessibility-policy>

## Share Your Story

The Harvard community has made this article openly available.

Please share how this access benefits you. [Submit a story](#)

WALKING POSTER PRESENTATION

Open Access

# Imaging surrogate of the ventricular arrhythmia substrate in a swine model of ventricular tachycardia: high resolution LGE vs. high resolution electro-anatomical mapping

Sébastien Roujol<sup>1\*</sup>, Tamer A Basha<sup>1</sup>, Cory M Tschabrunn<sup>1</sup>, Kraig V Kissinger<sup>1</sup>, Mark E Josephson<sup>1</sup>, Warren J Manning<sup>1,2</sup>, Elad Anter<sup>1</sup>, Reza Nezafat<sup>1</sup>

From 18th Annual SCMR Scientific Sessions  
Nice, France. 4-7 February 2015

## Background

Assessment of risk of sudden cardiac death (SCD) due to ventricular tachycardia and fibrillation (VT/VF) remains one of the most challenging areas in management of patients with chronic myocardial infarction. VT substrate consists of reentry circuits near or within a chronic scar. Invasive catheter mapping using electro-anatomical system (EAM) is the clinical standard for identification of reentry circuits. However, current EAM systems have low spatial resolution which may contribute to the high recurrence rate of VT after ablation (50% at 2 years). Better identification of reentry circuits can potentially improve SCD risk assessment and outcome of VT ablation procedures. Late gadolinium enhancement (LGE) CMR has the potential for non-invasive identification of reentry circuits, but LGE sequences suffer from low spatial resolution and partial voluming (1). In this study, we sought to utilize a high-resolution 3D LGE with 1 mm<sup>3</sup> isotropic spatial resolution to image the VT surrogate in swine model of VT.

## Methods

A 180 min balloon occlusion of the mid left anterior artery was induced in 6 Yorkshire swine. Each animal underwent an in-vivo CMR exam at 52±13 days after infarction. An electrophysiology study was subsequently performed using a novel proprietary high-resolution 64 electrode basket catheter (Rhythmia Medical) with programmed stimulation



to assess for VT inducibility and full LV mapping with ~7000 sampling points. Both in-vivo and ex-vivo CMR of the excised heart were performed using a 1.5 T Philips scanner. An injection of 0.2 mmol/kg of gadobenate dimeglumine was used during in-vivo CMR and again at 20 min before euthanasia. High resolution LGE was performed using a free breathing navigator-gated inversion recovery gradient echo (GRE) sequence (TR/TE/α=6.5/3.0ms/25°, FOV=270×270×112 mm<sup>3</sup>, voxel size=1×1×1 mm<sup>3</sup>, compressed sensing factor=4) and was reconstructed using a B<sub>1</sub>-weighted reconstruction based on LOST (2). Ex-vivo CMR included a T<sub>1</sub>-weighted acquisition with GRE imaging (TR/TE/α=15/7.1ms/25°, FOV=130×130×100 mm<sup>3</sup>, voxel size=0.5×0.5×0.5 mm<sup>3</sup>, 3 averages). In-vivo LGE data were analyzed offline. Endocardial and epicardial contours of the myocardium were manually delineated. The transmural intensity, averaged over the inner 1 mm layer of the myocardium, was normalized (where '0' corresponded to 3 standard deviations above the mean signal intensity of normal myocardium and '1' corresponds to 60% of the maximum signal intensity of the scar area) and projected on the endocardial shell.

## Results

Sustained reentrant VT was induced in all animals. In-vivo high resolution LGE shows good visual correspondence with ex-vivo CMR data (Figure 1). The inner 1 millimeter layer of the myocardium in LGE data was found to visually match with high resolution EAM (Figure 2).

<sup>1</sup>Department of Medicine (Cardiovascular Division), BIDMC / Harvard Medical School, Boston, MA, USA

Full list of author information is available at the end of the article

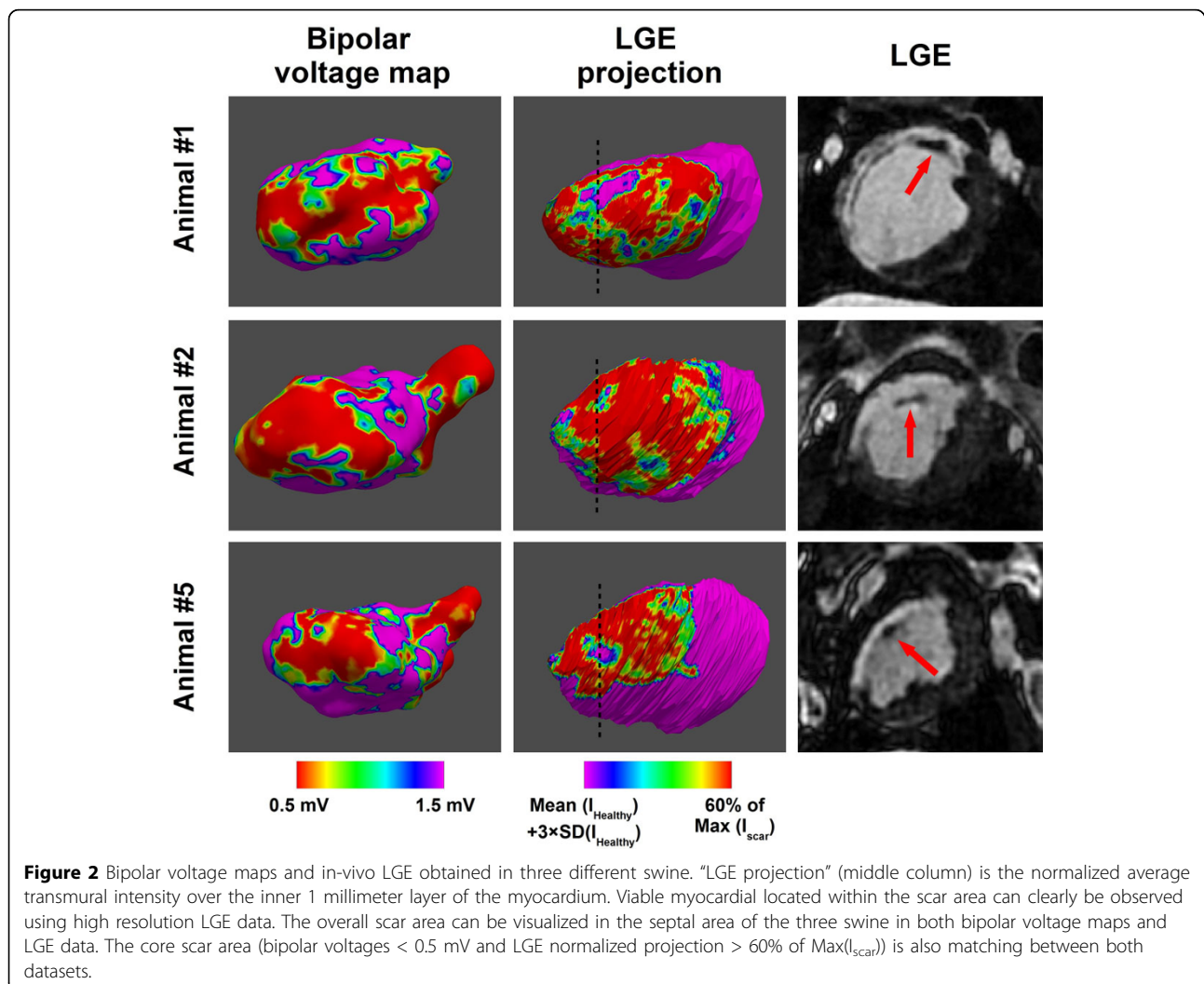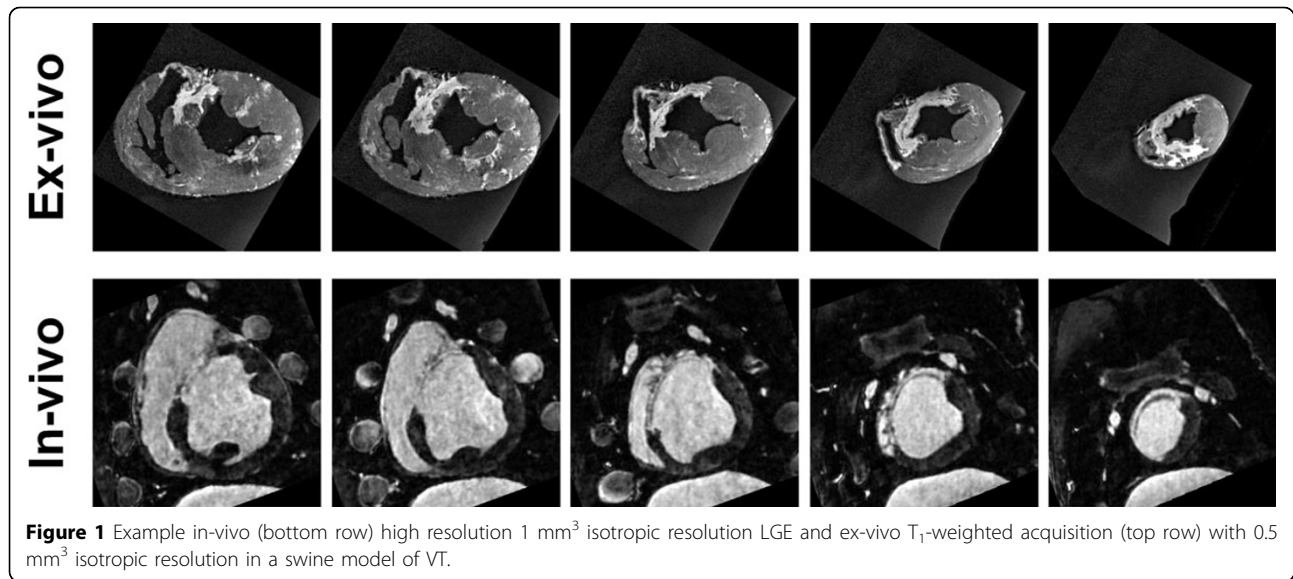

## Conclusions

High resolution LGE with 1 mm<sup>3</sup> isotropic resolution enables imaging of a surrogate of the VT substrate as defined by high-resolution EAM.

## Authors' details

<sup>1</sup>Department of Medicine (Cardiovascular Division), BIDMC / Harvard Medical School, Boston, MA, USA. <sup>2</sup>Radiology, BIDMC / Harvard Medical School, Boston, MA, USA.

Published: 3 February 2015

## References

1. Halperin: *JACC* 2011.
2. Akçakaya: *MRM* 2014.

doi:10.1186/1532-429X-17-S1-Q133

**Cite this article as:** Roujol et al.: Imaging surrogate of the ventricular arrhythmia substrate in a swine model of ventricular tachycardia: high resolution LGE vs. high resolution electro-anatomical mapping. *Journal of Cardiovascular Magnetic Resonance* 2015 **17**(Suppl 1):Q133.

**Submit your next manuscript to BioMed Central  
and take full advantage of:**

- Convenient online submission
- Thorough peer review
- No space constraints or color figure charges
- Immediate publication on acceptance
- Inclusion in PubMed, CAS, Scopus and Google Scholar
- Research which is freely available for redistribution

Submit your manuscript at  
[www.biomedcentral.com/submit](http://www.biomedcentral.com/submit)

