



Screening of DUB activity and specificity by MALDI-TOF mass spectrometry

Citation

Ritorto, M. S., R. Ewan, A. B. Perez-Oliva, A. Knebel, S. J. Buhrlage, M. Wightman, S. M. Kelly, et al. 2014. "Screening of DUB activity and specificity by MALDI-TOF mass spectrometry." *Nature communications* 5 (1): 4763. doi:10.1038/ncomms5763. <http://dx.doi.org/10.1038/ncomms5763>.

Published version

<https://doi.org/10.1038/ncomms5763>

Link

<http://nrs.harvard.edu/urn-3:HUL.InstRepos:14065303>

Terms of use

This article was downloaded from Harvard University's DASH repository, and is made available under the terms and conditions applicable to Other Posted Material (LAA), as set forth at

<https://harvardwiki.atlassian.net/wiki/external/NGY5NDE4ZjgzNTc5NDQzMGIzZWZhMGFIOWI2M2EwYTg>

Accessibility

<https://accessibility.huit.harvard.edu/digital-accessibility-policy>

Share Your Story

The Harvard community has made this article openly available. Please share how this access benefits you. [Submit a story](#)

Published in final edited form as:

Nat Commun. ; 5: 4763. doi:10.1038/ncomms5763.

Screening of DUB activity and specificity by MALDI-TOF mass spectrometry

Maria Stella Ritorto¹, Richard Ewan^{#1}, Ana B. Perez-Oliva^{#1}, Axel Knebel¹, Sara J. Buhrlage^{2,3}, Melanie Wightman¹, Sharon M. Kelly⁴, Nicola T. Wood¹, Satpal Virdee¹, Nathanael S. Gray^{2,3}, Nicholas A. Morrice⁵, Dario R. Alessi¹, and Matthias Trost^{1,*}

¹MRC Protein Phosphorylation and Ubiquitylation Unit, University of Dundee, Dundee, DD1 5EH, Scotland, UK

²Department of Cancer Biology, Dana-Farber Cancer Institute, Boston, MA 02115, U.S.A.

³Department of Biological Chemistry and Molecular Pharmacology, Harvard Medical School, 250 Longwood Avenue, SGM 628, Boston, MA 02115, U.S.A.

⁴Institute of Molecular Cell and Systems Biology, University of Glasgow, Glasgow, G12 8QQ, Scotland, UK

⁵The Beatson Institute for Cancer Research, Bearsden, Glasgow, G61 1BD, Scotland, UK

[#] These authors contributed equally to this work.

Abstract

Deubiquitylases (DUBs) are key regulators of the ubiquitin system which cleave ubiquitin moieties from proteins and polyubiquitin chains. Several DUBs have been implicated in various diseases and are attractive drug targets. We have developed a sensitive and fast assay to quantify *in vitro* DUB enzyme activity using matrix-assisted laser desorption/ionization time-of-flight (MALDI-TOF) mass spectrometry. Unlike other current assays, this method uses unmodified substrates, such as diubiquitin topoisomers. By analyzing 42 human DUBs against all diubiquitin topoisomers we provide an extensive characterization of DUB activity and specificity. Our results confirm the high specificity of many members of the OTU and JAMM DUB families and highlight that all USPs tested display low linkage selectivity. We also demonstrate that this assay can be deployed to assess the potency and specificity of DUB inhibitors by profiling 11 compounds against a panel of 32 DUBs.

Keywords

Deubiquitylase; ubiquitin; quantification; MALDI-TOF; mass spectrometry; inhibitor screening

*To whom correspondence should be addressed: **Matthias Trost**: Mailing address: MRC Protein Phosphorylation and Ubiquitylation Unit, College of Life Sciences, University of Dundee, Dow Street, Dundee, DD1 5EH, United Kingdom. Phone: +44 1382 386402. Fax: +44 1382 223778. m.trost@dundee.ac.uk.

Author contributions: MSR performed assays, all mass spectrometry experiments and analyzed the data, APO, RE, AK and SMK performed assays, RE and AK expressed and purified enzymes, SJB, NSG and SV provided compounds and expertise, MW and NTW performed cloning, MT, DRA, NAM and MSR designed research, MT, DRA and MSR wrote the manuscript with involvement of all other authors.

Conflict of interest statement: The authors declare no conflict of interest.

Post-translational modifications with ubiquitin control almost every process in cells. Ubiquitylation is facilitated by ubiquitin-activating (E1s), ubiquitin-conjugating (E2s), and ubiquitin ligase enzymes (E3s). Ubiquitin can be attached to substrate proteins as a single moiety or in the form of polymeric chains in which successive ubiquitin molecules are connected through specific isopeptide bonds. These bonds can be formed on any of the eight primary amines of the ubiquitin molecule (linear/N-terminus/M1, K6, K11, K27, K29, K33, K48 and K63) and thus can achieve a remarkable complexity, termed the ubiquitin code¹ in which the different chain topologies serve distinct signaling functions².

Ubiquitylation is reversible by specific cleavage through deubiquitylases (DUBs), of which about 90 have been identified in the human genome³. DUBs have been divided into five subclasses: Ubiquitin C-Terminal Hydrolases (UCHs), Ubiquitin-Specific Proteases (USPs), Machado-Joseph Disease Protein Domain Proteases (MJDs), Ovarian Tumor Proteases (OTUs), and JAB/MPN/Mov34 Metalloenzyme (JAMM) domain Proteases³⁻⁵. UCHs, USPs, OTUs and MJDs function as papain-like cysteine proteases, whereas JAMMs are zinc-dependent metalloproteases⁶. A sixth family of DUBs, Monocyte Chemotactic Protein Induced Proteases (MCPIP) has recently been proposed, but little is known about this family so far^{4, 6}.

DUBs have an essential role in ubiquitin homeostasis by catalyzing the editing and disassembly of polyubiquitin chains⁴. Furthermore, DUBs also perform signaling functions by the regulatory deubiquitylation of target proteins³ controlling proteasome-dependent protein degradation⁷, endocytosis⁸, DNA repair⁹ and kinase activation^{10, 11}. Not surprisingly, DUBs have been implicated in a number of diseases such as cancer¹²⁻¹⁷, inflammation^{10, 18}, neurodegeneration/Parkinson's disease¹⁹⁻²¹ and, due to their potentially drugable active sites, are considered attractive drug targets²².

Several chemical probes, such as Ub-vinyl methylester (VME), Ub-vinyl sulfone (VS)²³, branched and ubiquitin isopeptide activity-based probes²⁴ or diubiquitin activity probes²⁵ have been developed to explore the catalytic properties of DUBs. To screen for DUB inhibitors, current methods make use of non-physiological substrates including linear fusion of ubiquitin to a reporter protein such as phospholipase 2 or YFP in a Fluorescent Resonance Energy Transfer assay format^{26, 27}. Moreover, fusions of fluorogenic reporters such as Rhodamine²⁸ or 7-amino-4-methylcoumarin²⁹ to the C-terminal glycine of ubiquitin are also widely deployed. However, these substrates are not suitable for assessing the linkage specificity of DUBs. Furthermore, as these are artificial substrates that do not contain physiological isopeptide bonds, screening assays using these substrates could potentially identify compounds that might not inhibit the deubiquitylation of physiological substrates. To circumvent these issues it is possible to undertake DUB assays with more physiologically related diubiquitin molecules³⁰. However these assays are currently performed using low-throughput SDS-PAGE methodology and require relatively large amounts of enzymes (0.01-1 µg per assay) and substrates (typically up to 4 µg of substrate per assay)³¹.

Matrix-Assisted Laser Desorption/Ionization (MALDI) Time-of-Flight (TOF) mass spectrometry (MS)^{32, 33} has in the past been successfully applied to quantify low molecular

weight products of enzymes³⁴ or amyloid-beta peptides produced by gamma-secretase³⁵. Here, we present a novel screening method to assay DUB activity and specificity using unmodified diubiquitin isomer substrates. We employ quantitative MALDI-TOF MS using ¹⁵N-labeled ubiquitin and achieve high sensitivity, reproducibility and robustness. We analyze the specificity of 42 human DUBs and characterize the potency and selectivity of 11 DUB inhibitors against a panel of 32 DUBs. Our data represent an important resource for the scientific community and establish the applicability of the MALDI-TOF DUB assay in DUB inhibitor screening and selectivity assessment.

Results

MALDI-TOF DUB assay to assess DUB activity and specificity

We have developed a fast and sensitive assay to analyze *in vitro* activity and specificity of DUBs by MALDI-TOF mass spectrometry, termed the *MALDI-TOF DUB assay*. In this assay, we quantitate the amount of monoubiquitin generated by the *in vitro* cleavage of specific diubiquitin topoisomers by DUBs (Figure 1a). The DUB reaction consists of recombinant DUB (0.1–1,000 ng), diubiquitin (typically 125 ng/7,300 fmol) in 40 mM Tris-HCl pH 7.5, 5 mM dithiothreitol (DTT) and bovine serum albumin (BSA) carrier (0.25 µg) in a total volume of 5 µl. Reactions are undertaken for 1 hour at 30°C and terminated by addition of 1 µl of 10% (v/v) trifluoroacetic acid. Aliquots (2 µl) of each sample are spiked with 2 µl (1,000 fmol) of ¹⁵N-labeled ubiquitin (average mass 8666.55 Da), whose concentration was established by amino acid analysis, to serve as an internal standard for ubiquitin quantitation. A further 2 µl of 15.2 mg ml⁻¹ 2,5-dihydroxyacetophenone (DHAP) matrix and a 2 µl of 2% (v/v) trifluoroacetic acid are added and 0.5 µl of the resultant mixture are then spotted onto a 1,536 microtiter plate MALDI anchor target. The samples are analyzed by high mass accuracy MALDI-TOF MS in reflector positive ion mode on an UltrafleXtreme (Bruker Daltonics) mass spectrometer. The high resolution and mass accuracy of this MALDI-TOF mass spectrometer enabled baseline-resolution of isotopic patterns of ubiquitin and thus reliable quantitation of the area of the ubiquitin peak. Moreover, it permitted clear separation of the doubly-charged diubiquitin molecule (m/z 8556.64) and the singly-charged monoubiquitin (m/z 8565.76) (Figure 1b; Supplementary Figure 1). Next, we tested the linearity of our assay by analyzing standard curves over the ubiquitin concentration range of 10–10,000 nM (2–2,000 fmol on target) in the presence of 250 nM ¹⁵N-Ubiquitin (42 fmol on the target) and 874 nM diubiquitin (15 ng µl⁻¹; 146 fmol on the target) in three separately performed experiments on different days. Addition of ¹⁵N-ubiquitin and/or diubiquitin isomers, did not affect sensitivity with which ubiquitin could be detected and quantified (Supplementary Figure 2). Average correlation coefficient (r²) for the three curves was not less than 0.99 (Figure 1c) showing high linearity over a range of more than 500 (Supplementary Table 1). The mean intraday precision and interday accuracy for ubiquitin/¹⁵N-ubiquitin were 8% and 10% respectively, demonstrating the suitability of the assay as a screening tool. The lower limit of quantitation, defined as the lowest concentration that could be measured with a precision and accuracy better than 20%, was 10 nM (2 fmol on target) (Figure 1b) allowing for significantly reduced enzyme and substrate amounts compared to previously used low throughput methods that typically employed up to 4 µg of diubiquitin (234,000 fmol) per assay.

Determining DUB specificity

Utilizing the MALDI-TOF DUB assay, we systematically assessed the specificities of 42 recombinant human DUBs (Table 1) against all possible ubiquitin chain linkages. This represents almost 50% of the DUBs encoded in the human genome. For this, we determined the DUB activity at five different enzyme concentrations (from 0.02 to 200 ng μL^{-1}) against M1/linear, K6, K11, K27, K29, K33, K48 and K63-linked diubiquitin isomers, all at a final concentration of 1.46 μM in the assay. Altogether we performed more than 5,520 enzymatic reactions, providing the largest published resource for DUB specificity and activity (Figure 2). The results of this analysis highlighted a striking linkage specificity for five human DUBs (OTULIN – M1/linear, OTUB1 – K48, AMSH, AMSH-LP and BRCC3 – K63), which cleaved only one diubiquitin substrate even at very high concentrations of enzymes (Figure 2, **Group 1**), which is consistent with previous data analyzing these enzymes^{36–38} (Supplementary Table 2). Group 2 consisted of three DUBs that were highly specific to one linkage at only low concentrations (Cezanne-K11, OTUD1-K63 and A20-K48) and four DUBs (TRABID-K29/K33, VCPIPcat-K11/K48, OTUB2- and phosphorylated OTUD5-K48/K63) that displayed moderate selectivity hydrolyzing two ubiquitin linkages at low concentrations but were less selective at high concentrations (Figure 2, **Group 2**). Twenty DUBs, including all the active USP family members tested displayed little selectivity (Figure 2, **Group 3**), agreeing with previously reported findings³⁹. Four DUBs showed only very low activity (OTUD6A, OTU1, JOSD2 and ATXN3L) and six DUBs including OTU6B, JOSD1, and all the ubiquitin C-terminal hydrolases (UCH) were inactive in our assay (Figure 2, **Group 4**) (Supplementary Figure 3). In parallel, we also performed DUB fluorogenic Ubiquitin-110-Rhodamine-Glycine activity assays (Table 1) which are frequently used in the field^{28, 40, 41}. We calculated the specific activity of each DUB in this assay and grouped these into four categories (very low, low, moderate and high activity). Interestingly, when we compared the MALDI-TOF DUB assay data with fluorescent assay data (Table 1) we found 10 enzymes (USP9x, USP27x, USP36, CYLD, Otulin, OTUB1, OTUB2, AMSH, AMSH-LP and BRCC3) that were active only in the MALDI-TOF DUB assay. Four enzymes (USP10, USP28, A20 and VCPIP) displaying low activity in the fluorescence assay were significantly more active in the MALDI-TOF DUB assay. The majority (18 out of 42) of DUBs tested was active in the MALDI-TOF DUB assay and displayed moderate or high activity in the fluorescence assay. In contrast, seven DUBs including all members of the UCH family as well as OTUD6A, JOSD1 and JOSD2 were active in the Rhodamine assay but not in the MALDI-TOF DUB assay. We found that three DUBs tested (OTUD6B, OTU1, and ATXN3L) displayed very low or no activity, but nevertheless reacted with an active site directed probe (C-terminal propargylated ubiquitin)⁴² (Supplementary Figure 4).

Assessing potency and selectivity of DUB inhibitors

We next evaluated whether the MALDI-TOF DUB assay had potential to be deployed to assess potency and selectivity of DUB inhibitors. To undertake this, we set up a panel of 32 DUBs each assayed with the preferred diubiquitin isomers displaying the highest specific activity at the lowest concentration (Supplementary Figure 3). As proof-of-concept, we screened nine previously reported DUB inhibitors and inhibitor candidates (Compound 16⁴³,

L434078, WP1130⁴⁴, P22077⁴⁵, Febuxostat⁴⁶, SJB3-019A⁴⁷, PR-619⁴⁵, HBX 41,108⁴⁸, Pimozide⁴¹) as well as two E2/E3 ligase inhibitors that have potential to alkylate Cys residues (NSC 697923⁴⁹ and BAY 11-7082⁵⁰) at two concentrations against a group of 32 highly active DUBs from our assay. (Figure 3; Supplementary Figure 5; Supplementary Table 3). In addition, we performed IC₅₀ measurements for the DUBs that were most potently inhibited (Figure 4). For these studies, conditions were carefully optimized to ensure that assays were linear with respect to time (Supplementary Figure 6), and the diubiquitin substrate that displayed highest activity at the lowest DUB concentration was selected (Supplementary Figure 3).

Overall, none of the compounds tested displayed strong selectivity towards a single DUB and many were unselectively inhibiting most DUBs on the panel (Figure 3). For example, PR-619⁴⁵ (Supplementary Figure 5Q to R), is an ubiquitin/UbL isopeptidase inhibitor which has previously been reported to inhibit a range of cysteine protease DUBs⁴⁵. Consistent with this, we found that PR-619 inhibited 27 of the 32 tested DUBs with only members of the OTU and JAMM family being unaffected, likely because the latter are zinc-metalloproteases and do not possess a reactive catalytic Cys residue. Furthermore, our data indicates that SJB3-019A (Supplementary Figure 5O to P), which has previously been shown to inhibit USP1 in leukemic cells⁴⁷, inhibits USP8 more strongly (IC₅₀ 0.21 μ M) than USP1 (IC₅₀ 1.69 μ M) but in addition also significantly inhibited several other DUBs tested. Another reported USP1 inhibitor, Pimozide⁴¹, was also found to be non-selective, inhibiting many other DUBs with similar affinity to USP1 (Supplementary Figure 5E to F).

USP7 is one of the most targeted DUBs as phenotypes associated with USP7 silencing strongly suggest that small-molecule inhibitors of USP7 may have the potential for antiviral and anticancer therapies^{13, 51}. HBX 41,108, a cyano-indenopyrazine inhibitor of USP7 that has been shown to stabilize polyubiquitylated p53 at high concentrations in HEK293 cells⁴⁸, inhibited 25 of the 32 DUBs tested more than 70% at 5 μ M (Figure 3). This is consistent with another report suggesting that HBX 41,108 reacted with additional DUBs⁵². Only members of the JAMM family are not affected by HBX 41,108, again likely because they are zinc-metalloproteases and not cysteine proteases (Supplementary Fig. 5S to T).

Out of the 11 compounds profiled, BAY 11-7082 and NSC 697923 that both contain vinyl sulfone reactive groups (Supplementary Table 3) were found to inhibit USP7 with moderate higher potency than other DUBs tested. For example, BAY 11-7082 and NSC 697923 inhibited USP7 with an IC₅₀ of 0.19 μ M and 0.08 μ M, respectively. The next most potently inhibited DUB, i.e. USP21, was inhibited by BAY 11-7082 and NSC 697923 with an IC₅₀ of 0.96 μ M and 0.63 μ M, respectively. The potency of NSC 697923 for USP7 (0.08 μ M) was 70-fold higher than that of HBX 41,108 (5.97 μ M).

Discussion

We have developed a sensitive, reproducible and robust assay for the analysis of DUB *in vitro* activity and specificity. For this, we have made use of highly sensitive and fast MALDI-TOF mass spectrometry which is due to the use of 1,536-sample targets, suitable for robotic automation and thus high-throughput screening⁵³. We circumvented spot-to-spot

and shot-to-shot irreproducibility in MALDI ionization by using isotopically labeled ubiquitin as an internal standard as it guarantees identical extraction, crystallization and gas-phase behavior. Overall, this set-up allowed us to achieve very high precision, accuracy and linearity of measurements over concentrations of almost three orders of magnitude. The advantages compared to the commonly used assays with fluorogenic ubiquitin substrates are the use of substrates which are more physiological and the ability to analyze chain linkage specificity. Moreover, compared to current techniques using SDS-PAGE, our assay is considerably faster (2-4 hours for the acquisition of 1,536 data points) and more sensitive, thus requiring vastly reduced amounts of diubiquitin substrate. It should be noted that the assay is currently pipetted manually and due to addition of matrix and TFA, only 3.3% of the initial reaction mixture is utilized for the mass spectrometry analysis. Thus, after optimization, it should be feasible to scale down reaction amounts at least another 20-fold using nanoliter dispensing robotics representing a nearly 600-fold reduction in amounts of diubiquitin needed in current low throughput assays.

Our data has established that the MALDI-TOF DUB assay is a powerful approach to define the substrate specificity of DUBs. Using only 120 data points we have devised a strategy to characterize the activity of each DUB in triplicate (i.e. three different experiments) over five concentrations spanning 10,000-fold range against all eight diubiquitin chain linkages.

Only a few of the 42 expressed DUBs, and here particularly the members of the UCH family, were inactive in the MALDI-TOF DUB assay but showed high activity in the fluorogenic Ubiquitin-110-Rhodamine-Glycine assay (Table 1). This is consistent with previous work which has shown that UCH DUBs cleave ubiquitin moieties from protein substrates but do not hydrolyze diubiquitin^{45, 54}. Interestingly, we found that members of the JAMM family (AMSH, AMSH-LP, and BRCC3) displayed high activity in the MALDI-TOF DUB assay, exhibiting exquisite preference for K63 linkages, but were completely inactive in the fluorescence Ubiquitin-110-Rhodamine assay (Table 1). Therefore the MALDI-TOF DUB assay is the preferred technology to undertake future screening for specific inhibitors that target these metalloproteases. Among the other inactive DUBs are ATXN3L (MJD+ family DUB) that has been shown to preferentially cleave ubiquitin chains with more than four units³⁰ which is likely to explain why no activity was observed in the MALDI-TOF DUB and Ubiquitin-110-Rhodamine assays. Furthermore, OTUD6B was also previously shown to be inactive against ubiquitin dimers using a low-throughput assay³⁷. In our hands, full length OTU1 expressed in *E. coli* only displayed trace activity towards K11, K48 and K63 at the highest concentration tested (200 ng μl^{-1} ; 3.1 μM). In another study full length OTU1 assayed at 4 μM was shown to display low activity against K11, K27, K29, K33 and K48³⁷. A different group has suggested that OTU1 preferentially hydrolyzes longer polyubiquitin chains⁵⁵ which might explain the weak activity observed⁵⁵. Further work is also required to assess whether the other enzymes might require co-factors or post-translational modifications such as phosphorylation for optimal activity as reported for OTUD5⁵⁶. In the future, we intend to increase the coverage of the DUB family by including more enzymes. These proteins will be also expressed in either bacterial or insect cultures and if the full-length protein cannot be purified, a shorter construct encompassing the catalytic domain will be expressed.

Our data compares well to very recently published data of DUBs of the OTU family³⁷, confirming high specificity for many members of this family. Also other DUB families, such as the unspecific USPs as well as the specific JAMMs are in agreement with published data^{38, 39} (Supplementary Table 2). Yet, our data also suggests that the specificity of several DUBs depends on the concentration of the enzymes and the enzyme/substrate ratio. In general, highest selectivity is observed at low concentrations of DUBs. Cezanne for example, is very active and specific for K11 at the lowest concentration tested (i.e. 0.02 ng μl^{-1} ; ~ 0.2 nM). Similarly, OTUD1 is very selective for K63 at 0.2 ng μl^{-1} (~ 5 nM). However, at higher concentrations both Cezanne and OTUD1 lose their specificity. These observations highlight that these enzymes are not completely selective and possess the ability to weakly act on other topoisomers at higher substrate concentrations. Even several USPs, which are mostly unspecific in our assay, present some specificity at the lowest concentrations analyzed. This emphasizes the importance that specificity of DUBs should be tested over a wide range of enzyme concentrations, which has not generally been undertaken in previous analyses. The consistency of our data compared with previous work on DUB activity and selectivity highlights the reliability of the MALDI-TOF DUB assay technology.

None of the DUBs tested initially displayed significant activity against K27-linked diubiquitin isomers that were purchased commercially. We confirmed by mass spectrometry that the commercial K27 diubiquitin molecule was indeed correctly linked and was present in equimolar amounts compared to the other diubiquitin isomers. We determined by parallel reaction monitoring against the linkage peptides that this diubiquitin contained small amounts ($\sim 7\%$) of K11 and ($\sim 8\%$) K63 diubiquitin chains which were not affecting overall results though (Supplementary Figures 10 and 11A). Additionally, we performed circular dichroism (CD) to rule out misfolding of the isomer (Supplementary Figure 11B). However, to ensure that the inactivity of all DUBs against this linkage was not due to quality issues of the commercially produced K27 diubiquitin, we compared three differently sourced versions of K27 diubiquitin (Boston Biochem, UbiQ and in-house chemically engineered K27 diubiquitin^{57, 58}) against four DUBs that had shown activity against K27 in previous publications^{25, 39}. While USP7, USP8 and USP25 could not hydrolyze any of the three differently sourced K27 diubiquitins in our hands, we observed that USP16 displayed strong activity against this linkage for in-house K27 diubiquitin but not for the commercially sourced K27-chains (Supplementary Figure 11C). It should be noted that previous work³⁹ also concluded that USP16 displayed highest activity towards K27-linked diubiquitin isomers compared to the other DUBs tested, suggesting that this enzyme might indeed cleave this chain type *in vitro*. It is not clear why there is this discrepancy in activity of USP16 towards in-house chemically engineered (GOPAL methodology) and commercially available K27 diubiquitin. We therefore re-tested the whole DUB panel against the in-house K27 diubiquitin and except for USP16, no other DUB showed any significant activity, agreeing with the data obtained with commercially sourced K27 diubiquitin. This suggests that extra caution is required when K27 diubiquitin is used for DUB assays.

There is huge interest in developing chemical probes that target specific components of the ubiquitylation system⁵⁹. We have shown that the MALDI-TOF DUB assay can be readily used to determine inhibition rates and the IC_{50} of small molecule inhibitors of DUBs. The

MALDI-TOF DUB assay also enables the facile profiling of inhibitors against numerous DUBs acting on a more physiological substrate than fluorescent ubiquitin conjugates that have been previously employed for this purpose²⁸. Moreover, one will be able to employ this assay using other physiological substrates, such as ubiquitylated proteins. As proof-of-concept, we have deployed a panel of 32 active DUBs to profile 11 available DUB inhibitors and inhibitor candidates. Our work confirms previous work that PR-619 is a general DUB inhibitor that potently suppresses the activity of almost all cysteine protease DUBs⁴⁵. Similarly HBX 41,108, proposed as an USP7 inhibitor⁴⁸, inhibited almost all DUBs in our assay better than USP7. Out of the compounds analyzed, BAY 11-7082 and NSC 697923 displayed the highest selectivity, inhibiting USP7 with 5-8 fold higher potency than the second most sensitive DUB on our panel (i.e. USP21). BAY 11-7082 inhibits NF κ B signaling⁶⁰ and was recently shown to inhibit the majority of E2 and E3 ligases tested by reacting covalently with the catalytic cysteine residues⁵⁰. Moreover, BAY 11-7082 also inhibits several tyrosine phosphatases by reacting with catalytic Cys residue of these enzymes⁶¹. NSC 697923 was originally shown to inhibit the E2 ligase Ubc13-Uev1A⁴⁹. These data suggest that BAY 11-7082 and NSC 697923 are likely to inhibit a broad range of enzymes possessing catalytic Cys residues. Nevertheless, the moderate specificity of these compounds for USP7 indicates that it might be possible to further engineer these compounds towards more selective probes. Furthermore, we could show that SJB3-019A and Pimozide, two proposed USP1 inhibitors^{41, 47}, inhibit USP8 at a ten-fold lower concentration or showed poorly selectivity, respectively. Overall this data reveals the importance of undertaking extensive profiling of specificity of DUB inhibitors as it is essential to ensure the selectivity of these compounds for *in vivo* applications.

In order to screen larger numbers of small molecule inhibitors, we propose a future screening strategy (Supplementary Figure 7) that can be summarized in three steps: Screen 1, screening a large number of small molecules against a single DUB to identify lead candidates; Screen 2, inhibitor specificity determination of these lead candidates at a single concentration against a panel of a large number of DUBs; Screen 3, determination of IC₅₀ for best DUB/inhibitor pairs. We believe that such strategies will be useful for discovery of specific inhibitors of DUBs which have the potential to become important future drug targets²².

In conclusion, we present here a novel screening method to assay DUB activity and specificity with high sensitivity, reproducibility and reliability, which is able to carry out precise quantitative measurements at a rate of ~6-9 s per sample spot. Using physiological substrates of DUBs allowed us to determine specificity of 42 human DUBs among which several showed high specificity for one single chain type. This data allowed us to generate a simple array of preferred chain types and lowest concentrations of activity for each DUB which will serve as a sensitive and fast tool for screening for DUB inhibitors.

Methods

Materials

Ubiquitin monomer, Bovine Serum Albumin (BSA), Tris and DTT were purchased from Sigma Aldrich. Diubiquitin topoisomers (M1, K6, K11, K27, K29, K33, K48 and K63-

linked) were purchased from Boston Biochem (Boston, MA), additional K27 diubiquitin was produced in-house^{57, 58}, whereas all MALDI-TOF MS materials (targets, matrix and protein calibration mixture) were purchased from Bruker Daltonics (Bremen, Germany). PR-619 and P22077 (Calbiochem/Merck, Darmstadt, Germany) as well as HBX 41,108, Pimozide and Degrasyn/WP1130 (Tocris Bioscience, Bristol, UK) and L434078 (Sigma-Aldrich, St. Louis, MO, USA) were purchased commercially. Febuxostat, SJB3-019A, Compound 16, NSC 697923 and BAY 11-7082 were synthesized (Supplementary Table 3).

Expression of DUB enzymes

For bacterial expression, full length or catalytic core domains for various DUB enzymes were cloned into either pET-24, pET-28 (Novagen) or pGEX6P (GE-Healthcare) vectors to express either N-terminally tagged 6×His or GST tagged proteins. Recombinant proteins were expressed in *E. coli* Bl21 (DE3) cells which were lysed by sonication lysis buffer (50 mM Tris-HCl pH 7.5, 250 mM NaCl, 0.5 mM EGTA, 0.5 mM EDTA, 0.5% Triton-X100, 1 mM DTT, 1mM Pefabloc, 10 µg ml⁻¹ Leupeptin), then centrifuged to remove insoluble material. For protein purification, supernatants were subjected to affinity chromatography using either Ni²⁺-NTA-Sepharose (GE Healthcare) or Glutathione-Sepharose (Expediton) resin. For insect cell protein expression, appropriate cDNAs were cloned into the pFastBac vector, baculoviruses were generated to encode various Dac-tagged⁶² DUB enzymes. Sf21 cells were typically infected with P1 virus stocks and harvested 48 hours later. Cells were lysed in Dac lysis buffer (40 mM Tris pH 7.5, 0.2% Triton-X 100, 0.5 mM EGTA, 0.1 mM EDTA, 1mM DTT) supplemented with 1 mM Pefabloc and Leupeptin at 20 µg ml⁻¹ then centrifuged to remove insoluble material. For protein purification, supernatants were subjected to affinity chromatography using ampicillin-Sepharose resin for 45 min at ambient temperature. The DUB enzymes were either eluted from the ampicillin-Sepharose by incubating 4 × for 15min with 1 resin volume of 50 mM Tris-HCl pH 7.5, 5% v/v glycerol, 100 mM NaCl, 10 mM ampicillin, 1 mM DTT, 0.03% (w/v) Brij-35 or recovered by digesting the DUB off the Dac-tag using TEV-protease in 50 mM Tris-HCl pH 7.5, 100 mM NaCl, 0.03% (w/v) Brij 35 (Supplementary Figure 8).

¹⁵N-ubiquitin expression and purification

Untagged full length human ubiquitin was cloned into the pET24 vector and expressed in *E. coli* Bl21 (DE3) cells. Cells expressing ¹⁵N-ubiquitin were grown in M9 minimal media supplemented with ¹⁵N ISOGRO (Sigma) and Ammonium-¹⁵N chloride (Sigma) according to manufacturer's instructions. Cells were sedimented, resuspended in H₂O and lysed by freeze-thawing, then centrifuged to remove insoluble material. Bacterial proteins were precipitated by dropping the pH to 4.5 with diluted HClO₄, then sedimented by centrifugation. The supernatant containing the ubiquitin was adjusted to 20 mM ammonium-acetate pH 4.5 and applied to a Source 15 S HR10/10 column (GE Healthcare), which was developed with a gradient of 0-1 M NaCl. The ¹⁵N-ubiquitin eluting at a conductivity of 18 mS cm⁻¹ was concentrated and subjected to chromatography on a Superdex75 XK16/60 column (Amersham). Final ¹⁵N-ubiquitin fractions were pooled and concentrated to 35 mg ml⁻¹. MALDI-TOF MS analysis revealed a 97% incorporation of ¹⁵N (Supplementary Figure 9).

Preparation of K27-linked diubiquitin

K27-linked diubiquitin was prepared as described⁵⁷ with the following exceptions. Acceptor ubiquitin was expressed from a new plasmid (pCDF-pylT-UbTAG27-His6) carrying an amber stop codon at residue position 27 of ubiquitin. Crude Cbz-deprotected ubiquitin species were then dissolved in denaturing buffer (200 mM Na₂HPO₄ pH 7.5, 6 M guanidinium chloride) at a final concentration of ~2mg ml⁻¹. Deprotected diubiquitin was then purified from residual mono ubiquitin by semi-preparative reversed-phase HPLC using a Dionex Ultimate 3000 system. A flow rate of 10 ml min⁻¹ and a gradient of 20 % solvent A to 50 % solvent B over 40 minutes were used with a Thermo Biobasic C4 (250 mm × 21.2 mm) column (solvent A = 0.1 % trifluoroacetic acid in H₂O; solvent B = 0.1 % trifluoroacetic acid in acetonitrile). Fractions corresponding to K27-linked diubiquitin were determined by LC-MS and were then pooled and freeze dried. Freeze dried K27-linked diubiquitin was dissolved in denaturing buffer to a final concentration of 2 mg ml⁻¹ and folded by overnight dialysis against PBS.

In vitro DUB assays and inhibitor screening

To monitor DUB activity *in vitro*, we tested a panel of 42 human DUBs at different concentrations (0.02-0.2-2-20-200 ng µl⁻¹) against all diubiquitin topoisomers (M1, K6, K11, K27, K29, K33, K48 and K63-linked chains) as substrates at a fixed concentration (25 ng µl⁻¹, 1.46 µM). Both enzymes and substrates were freshly prepared in the reaction buffer (40 mM Tris-HCl, pH 7.6, 5 mM DTT, 0.005% BSA) for each run. The enzymes were pre-incubated in the reaction buffer for 10 min at 30°C; afterwards, the diubiquitin isomers were added and the reaction mixture incubated for 60 min at 30°C. The reaction was stopped by adding TFA to a final concentration of 2% (v/v). Possible background due to contamination of the diubiquitin with ubiquitin monomers was measured in a reaction buffer in which the enzyme was excluded and ubiquitin intensities normalized accordingly (Supplementary Figure 10B-2). Dimer purity was controlled by SDS-PAGE, PRM and MALDI-TOF MS/MS (Supplementary Figures 10-11).

For the small molecule inhibitor studies, we tested the MALDI-TOF DUB methodology by screening 11 DUB inhibitors or inhibitor candidates. To assess linearity for inhibitor experiments, time-dependent inhibition experiments were performed by adding increasing concentrations of the compound, from 0 to 90 µM, to the reaction buffer containing USP7 (1 ng µl⁻¹) and incubated for 30 min at 30°C. We demonstrated the linearity of response of USP7 (1 ng µl⁻¹) *versus* K11-linked diubiquitin (1.46 µM) to increasing concentrations of PR-619 (0-40 µM) over time (0-20 min) (Supplementary Figure 6). Some of these inhibitors could potentially react with reducing agents present in the assay such as dithiothreitol (DTT) or tris(2-carboxyethyl)phosphine (TCEP). We therefore verified using MALDI-TOF mass spectrometry whether the inhibitor compounds reacted with either 1 mM DTT or 0.5 mM TCEP when incubated at 30°C for 1 hour under the DUB assay conditions employed. This revealed that only WP1130 significantly reacted with both DTT and TCEP under these conditions (data not shown). For all inhibitor experiments, except PR-619 and HBX 41,108, no DTT but the remaining trace levels from the protein expression was added to the reaction buffer.

Following the scheme in Supplementary Figure 7, we determined the activity of the inhibitors in two steps. First of all, we screened the inhibitors in duplicates against 32 DUBs. The enzymes were pre-incubated with either 1 or 3/5/10 μM of inhibitor for 35 min at 30°C. Next, substrates were incubated with the enzyme plus inhibitor mixture for 60 min at 30°C. Secondly, we determined the inhibitors' IC_{50} on a subset of selected DUBs. For calculation of IC_{50} , data were fitted in SigmaPlot (v. 12.5 Build 12.5.0.38) using the four parameter logistic equation: $y = \text{min} + (\text{max} - \text{min}) / (1 + (x / \text{IC}_{50})^{\text{Hillslope}})$. Values of IC_{50} for all compounds are summarized in Supplementary Table 4.

Analysis by MALDI-TOF MS

Sample preparation—Acidified samples of the DUB assays were mixed with 0.5 μM ^{15}N -ubiquitin and then with one part of 2 % (v/v) TFA and one part of DHAP matrix solution (7.6 mg of DHAP in 375 μl ethanol and 125 μl of an aqueous 12 mg ml^{-1} diammonium hydrogen citrate). 0.5 μl of these solutions were spotted in replicates onto an MTP AnchorChip 1,536 TF (600 μm anchor, Bruker Daltonics).

Data acquisition—A high resolution MALDI-TOF MS instrument (UltrafleXtreme, Bruker Daltonics) with Compass 1.3 control and processing software was used. Samples were run in automatic mode (AutoXecute, Bruker Daltonics) allowing 6-9 seconds per spot, using the 1,536 spots AnchorChip. Ionization was achieved by a 1 kHz smartbeam-II solid state laser with a fixed initial laser power of 60% (laser attenuator offset 68%, range 30%) and detected by the FlashDetector at detector gain of 10 $\times$. Reflector mode was used with optimized voltages for reflector-1 (26.61 kV) and reflector-2 (13.39 kV), ion sources (IonSource-1: 24.86 kV, IonSource-2: 22.71 kV) and pulsed ion extraction (320ns). Sampling rate was 0.25 ns equivalent to a 4 GS/s digitization rate. An amount of 2100 (3 $\times$ 700) shots were summed up in “random walk” and with “large” smartbeam laser focus. Spectra were accumulated by FlexControl software (v. 3.3 Build 108), processed using FlexAnalysis software (v. 3.3, Build 80) and the Sophisticated Numerical Annotation Procedure (“SNAP”) peak detection algorithm, setting the signal to noise threshold at 250. Prior to calibration, the spectra were processed using smoothing (Savitzky-Golay algorithm) and baseline subtraction (“TopHat”) for reproducible peak annotation on non-resolved isotope distributions: 1 cycle, 0.2 m/z for the width. For external interactive calibration in quadratic mode, the “Protein Calibration Standard 1” (Bruker) was used with ubiquitin ($[\text{M} + \text{H}]^+ \text{avg} = 8565.76$), myoglobin ($[\text{M} + 2\text{H}]^{2+} \text{avg} = 8476.66$) and cytochrome c ($[\text{M} + 2\text{H}]^{2+} \text{avg} = 6181.05$, $[\text{M} + \text{H}]^+ = 12360.97$) ions as average m/z values. Internal calibration was performed using the ubiquitin peak ($[\text{M} + \text{H}]^+ \text{avg} = 8565.76$). An example of a full scan MALDI spectrum is depicted in Supplementary Figure 12.

Data analysis—A modified method for data acquisition was developed for FlexAnalysis Software, using the SNAP algorithm. This algorithm is calculating an isotopic distribution which best fits the pattern of an isotopically resolved protein signal. For area calculation, the complete isotopic distribution was taken into account. Data output was exported as a .csv file using FlexAnalysis Batch Process (Compass 1.3) and further processed in Microsoft Excel, where plotting of graphs, calculation of standard deviation (SD), and coefficient of variation CV (%) were performed. The measurement of DUB activity (% of diubiquitin

isomer consumed, x) from relative isotopic distribution summed area ratios was performed according to equation (1).

$$x = \left(\left(\frac{\text{Area}_{Ubi}}{\text{Area}_{^{15}NUbi}} * [^{15}NUbi] \right) / 2 \right) / [DiUbi] * 100 \quad \text{equation (1)}$$

Quality control of diubiquitin isomers by MALDI-TOF and PRM

Two μg of each diubiquitin chain was resolved by SDS-PAGE (4-20%, Tris-Glycine, Novex, Life Technologies), stained by InstantBlue coomassie stain (Expedeon). Bands were excised and digested with trypsin (Supplementary Figure 10A). Gel pieces were washed subsequently with water, 50% methanol/water, 0.1M NH_4HCO_3 , then shrunk in acetonitrile and digested with trypsin (Pierce). The digestion was performed overnight at 37°C and peptides were extracted further in 50% acetonitrile/2.5% TFA. Digests were dried and reconstituted in 0.1% TFA/water to 15 $\mu\text{g ml}^{-1}$. A total of 30 ng of each digest was injected onto a 15 cm $\times$ 75 μm (I.D.) EasySpray column (Thermo Fisher Scientific) and analyzed on a LTQ-Velos Pro ion trap (Thermo Fisher Scientific) using a parallel reaction monitoring program (PRM)⁶³ specific for each diubiquitin chain linkage (Gly-Gly) peptide (Supplementary Figures 10C and 11A). Lists of precursor masses and fragment transitions are reported in Supplementary Table 5. Data was acquired in a data-independent mode with one full scan followed by 10 MS^2 scans with the masses of the different linkages. MS^2 occurred even if precursor mass was not detected in MS^1 scan. Extracted ion chromatograms of the MS^2 spectra for each diubiquitin chain peptide was performed by summing the ion current of the three or four most dominant daughter ions.

The linkage peptide of K29-diubiquitin does not bind to the trap column (Acclaim PepMap100 C18, 5 μm , Thermo Fisher Scientific) of the LC-MS system under normal conditions. We therefore analyzed the purity of this diubiquitin by MALDI-TOF MS (Supplementary Figure 10D). Fifteen $\text{ng } \mu\text{l}^{-1}$ of digested diubiquitin was mixed with 10 mg ml^{-1} *alpha*-cyano-4-hydroxycinnamic acid (1:1) and spotted onto a 384 AnchorChip target (Bruker Daltonics). For MS^2 , LIFT technology was performed and the data were processed by Mascot server through BioTools (Bruker Daltonics).

Ubiquitin-Rhodamine assay

Ubiquitin-Rhodamine-110-Glycine was prepared in-house²⁸. 0.5 μM Ubiquitin-Rhodamine-110-Glycine in 40 mM Tris-HCl buffer, pH 7.6, 5 mM DTT and 0.05 mg ml^{-1} BSA were incubated with 0.05-5 $\text{ng } \mu\text{l}^{-1}$ of each DUB for 60 min at 30°C. Samples were prepared in triplicates and analyzed in 96-well plates using a Perkin Elmer Envision 2104 multi label reader at Excitation/Emission 485/535 nm^{28} .

Supplementary Material

Refer to Web version on PubMed Central for supplementary material.

Acknowledgments

This work was funded by Medical Research Council UK and the pharmaceutical companies supporting the Division of Signal Transduction Therapy (DSTT) (AstraZeneca, Boehringer-Ingelheim, GlaxoSmithKline, Janssen Pharmaceutica, Merck KGaA and Pfizer). We would like to thank the DNA cloning, Protein Production, DNA sequencing and Mass Spectrometry teams of the MRC Protein Phosphorylation and Ubiquitylation Unit for their support, Yogesh Kulathu for expression of vOTU and suggestions and Natalia Shpiro for the synthesis of compounds. We also thank Bruker Daltonics and particularly Anja Resemann and Rainer Paape, as well as Rod Watson and Julia Smith for providing scripts and technical support.

ABBREVIATIONS

AVG	average
BSA	bovine serum albumin
DHAC	diammonium hydrogen citrate
DHAP	2,5 dihydroxyacetophenone
DUB	Deubiquitylase
HTS	high throughput screening
DTT	dithiothreitol
JAMM	JAB/MPN/Mov34 Metalloenzyme
MALDI	matrix-assisted laser desorption/ionization
MCPIP	Monocyte Chemotactic Protein Induced Proteases
MJD	Machado-Joseph Disease Protein Domain Protease
MS	mass spectrometry
MTP	microtiter plate
OTU	Ovarian Tumor Protease
TOF	time-of-flight
UCH	Ubiquitin C-Terminal Hydrolase
USP	Ubiquitin Specific Protease

References


1. Komander D, Rape M. The ubiquitin code. *Annu. Rev. Biochem.* 2012; 81:203–229. [PubMed: 22524316]
2. Ikeda F, Dikic I. Atypical ubiquitin chains: new molecular signals. ‘Protein Modifications: Beyond the Usual Suspects’ review series. *EMBO Rep.* 2008; 9:536–542. [PubMed: 18516089]
3. Nijman SM, et al. A genomic and functional inventory of deubiquitinating enzymes. *Cell.* 2005; 123:773–786. [PubMed: 16325574]
4. Komander D, Clague MJ, Urbe S. Breaking the chains: structure and function of the deubiquitinases. *Nat. Rev. Mol. Cell Biol.* 2009; 10:550–563. [PubMed: 19626045]
5. Liang J, et al. MCP-induced protein 1 deubiquitinates TRAF proteins and negatively regulates JNK and NF-kappaB signaling. *J. Exp. Med.* 2010; 207:2959–2973. [PubMed: 21115689]
6. Reyes-Turcu FE, Ventii KH, Wilkinson KD. Regulation and cellular roles of ubiquitin-specific deubiquitinating enzymes. *Annu. Rev. Biochem.* 2009; 78:363–397. [PubMed: 19489724]

7. Lee MJ, Lee BH, Hanna J, King RW, Finley D. Trimming of ubiquitin chains by proteasome-associated deubiquitinating enzymes. *Mol. Cell. Proteomics*. 2011; 10:R110 003871. [PubMed: 20823120]
8. Clague MJ, Urbe S. Endocytosis: the DUB version. *Trends Cell Biol*. 2006; 16:551–559. [PubMed: 16996268]
9. Nijman SM, et al. The deubiquitinating enzyme USP1 regulates the Fanconi anemia pathway. *Mol. Cell*. 2005; 17:331–339. [PubMed: 15694335]
10. Kovalenko A, et al. The tumour suppressor CYLD negatively regulates NF-kappaB signalling by deubiquitination. *Nature*. 2003; 424:801–805. [PubMed: 12917691]
11. Al-Hakim AK, et al. Control of AMPK-related kinases by USP9X and atypical Lys(29)/Lys(33)-linked polyubiquitin chains. *Biochem. J*. 2008; 411:249–260. [PubMed: 18254724]
12. Zhang Y, et al. USP44 regulates centrosome positioning to prevent aneuploidy and suppress tumorigenesis. *J. Clin. Invest*. 2012; 122:4362–4374. [PubMed: 23187126]
13. Li M, et al. Deubiquitination of p53 by HAUSP is an important pathway for p53 stabilization. *Nature*. 2002; 416:648–653. [PubMed: 11923872]
14. Graner E, et al. The isopeptidase USP2a regulates the stability of fatty acid synthase in prostate cancer. *Cancer Cell*. 2004; 5:253–261. [PubMed: 15050917]
15. Oliveira AM, et al. USP6 (Tre2) fusion oncogenes in aneurysmal bone cyst. *Cancer Res*. 2004; 64:1920–1923. [PubMed: 15026324]
16. Shah SP, et al. Mutational evolution in a lobular breast tumour profiled at single nucleotide resolution. *Nature*. 2009; 461:809–813. [PubMed: 19812674]
17. van der Horst A, et al. FOXO4 transcriptional activity is regulated by monoubiquitination and USP7/HAUSP. *Nat. Cell Biol*. 2006; 8:1064–1073. [PubMed: 16964248]
18. Harhaj EW, Dixit VM. Deubiquitinases in the regulation of NF-kappaB signaling. *Cell Res*. 2011; 21:22–39. [PubMed: 21119682]
19. Liu Y, Fallon L, Lashuel HA, Liu Z, Lansbury PT Jr. The UCH-L1 gene encodes two opposing enzymatic activities that affect alpha-synuclein degradation and Parkinson's disease susceptibility. *Cell*. 2002; 111:209–218. [PubMed: 12408865]
20. Saigoh K, et al. Intragenic deletion in the gene encoding ubiquitin carboxy-terminal hydrolase in gad mice. *Nat. Genet*. 1999; 23:47–51. [PubMed: 10471497]
21. Kondapalli C, et al. PINK1 is activated by mitochondrial membrane potential depolarization and stimulates Parkin E3 ligase activity by phosphorylating Serine 65. *Open Biol*. 2012; 2:120080. [PubMed: 22724072]
22. Cohen P, Tcherpakov M. Will the ubiquitin system furnish as many drug targets as protein kinases? *Cell*. 2010; 143:686–693. [PubMed: 21111230]
23. Hemelaar J, et al. Specific and covalent targeting of conjugating and deconjugating enzymes of ubiquitin-like proteins. *Mol. Cell. Biol*. 2004; 24:84–95. [PubMed: 14673145]
24. Iphof A, et al. Profiling ubiquitin linkage specificities of deubiquitinating enzymes with branched ubiquitin isopeptide probes. *Chembiochem*. 2012; 13:1416–1420. [PubMed: 22689415]
25. McGouran JF, Gaertner SR, Altun M, Kramer HB, Kessler BM. Deubiquitinating enzyme specificity for ubiquitin chain topology profiled by di-ubiquitin activity probes. *Chem. Biol*. 2013; 20:1447–1455. [PubMed: 24290882]
26. Nicholson B, et al. Characterization of ubiquitin and ubiquitin-like-protein isopeptidase activities. *Protein Sci*. 2008; 17:1035–1043. [PubMed: 18424514]
27. Horton RA, Strachan EA, Vogel KW, Riddle SM. A substrate for deubiquitinating enzymes based on time-resolved fluorescence resonance energy transfer between terbium and yellow fluorescent protein. *Anal. Biochem*. 2007; 360:138–143. [PubMed: 17118327]
28. Hassiepen U, et al. A sensitive fluorescence intensity assay for deubiquitinating proteases using ubiquitin-rhodamine110-glycine as substrate. *Anal. Biochem*. 2007; 371:201–207. [PubMed: 17869210]
29. Dang LC, Melandri FD, Stein RL. Kinetic and mechanistic studies on the hydrolysis of ubiquitin C-terminal 7-amido-4-methylcoumarin by deubiquitinating enzymes. *Biochemistry*. 1998; 37:1868–1879. [PubMed: 9485312]


30. Nicastro G, et al. Understanding the role of the Josephin domain in the PolyUb binding and cleavage properties of ataxin-3. *PLoS ONE*. 2010; 5:e12430. [PubMed: 20865150]
31. Hospenthal MK, Freund SM, Komander D. Assembly, analysis and architecture of atypical ubiquitin chains. *Nat. Struct. Mol. Biol.* 2013
32. Karas M, Hillenkamp F. Laser desorption ionization of proteins with molecular masses exceeding 10,000 daltons. *Anal. Chem.* 1988; 60:2299–2301. [PubMed: 3239801]
33. Szajli E, Feher T, Medzihradszky KF. Investigating the quantitative nature of MALDI-TOF MS. *Molecular and cellular proteomics*. 2008; 7:2410–2418. [PubMed: 18653768]
34. Bungert D, Heinze E, Tholey A. Quantitative matrix-assisted laser desorption/ionization mass spectrometry for the determination of enzyme activities. *Anal. Biochem.* 2004; 326:167–175. [PubMed: 15003557]
35. Rufenacht P, Guntert A, Bohrmann B, Ducret A, Dobeli H. Quantification of the A beta peptide in Alzheimer's plaques by laser dissection microscopy combined with mass spectrometry. *Journal of mass spectrometry: JMS*. 2005; 40:193–201. [PubMed: 15706631]
36. Keusekotten K, et al. OTULIN antagonizes LUBAC signaling by specifically hydrolyzing Met1-linked polyubiquitin. *Cell*. 2013; 153:1312–1326. [PubMed: 23746843]
37. Mevissen TE, et al. OTU deubiquitinases reveal mechanisms of linkage specificity and enable ubiquitin chain restriction analysis. *Cell*. 2013; 154:169–184. [PubMed: 23827681]
38. McCullough J, Clague MJ, Urbe S. AMSH is an endosome-associated ubiquitin isopeptidase. *The Journal of cell biology*. 2004; 166:487–492. [PubMed: 15314065]
39. Faesen AC, et al. The differential modulation of USP activity by internal regulatory domains, interactors and eight ubiquitin chain types. *Chem. Biol.* 2011; 18:1550–1561. [PubMed: 22195557]
40. Orcutt SJ, Wu J, Eddins MJ, Leach CA, Strickler JE. Bioluminescence assay platform for selective and sensitive detection of Ub/Ubl proteases. *Biochim. Biophys. Acta*. 2012; 1823:2079–2086. [PubMed: 22705352]
41. Chen J, et al. Selective and cell-active inhibitors of the USP1/ UAF1 deubiquitinase complex reverse cisplatin resistance in non-small cell lung cancer cells. *Chem. Biol.* 2011; 18:1390–1400. [PubMed: 22118673]
42. Ekkebus R, et al. On terminal alkynes that can react with active-site cysteine nucleophiles in proteases. *J. Am. Chem. Soc.* 2013; 135:2867–2870. [PubMed: 23387960]
43. Colombo M, et al. Synthesis and biological evaluation of 9-oxo-9H-indeno[1,2-b]pyrazine-2,3-dicarbonitrile analogues as potential inhibitors of deubiquitinating enzymes. *Chemmedchem*. 2010; 5:552–558. [PubMed: 20186914]
44. Pham LV, et al. Degrasyn potentiates the antitumor effects of bortezomib in mantle cell lymphoma cells in vitro and in vivo: therapeutic implications. *Mol. Cancer Ther.* 2010; 9:2026–2036. [PubMed: 20606045]
45. Altun M, et al. Activity-based chemical proteomics accelerates inhibitor development for deubiquitylating enzymes. *Chem. Biol.* 2011; 18:1401–1412. [PubMed: 22118674]
46. Takano Y, et al. Selectivity of febuxostat, a novel non-purine inhibitor of xanthine oxidase/ xanthine dehydrogenase. *Life Sci.* 2005; 76:1835–1847. [PubMed: 15698861]
47. Mistry H, et al. Small-molecule inhibitors of USP1 target ID1 degradation in leukemic cells. *Mol. Cancer Ther.* 2013; 12:2651–2662. [PubMed: 24130053]
48. Colland F, et al. Small-molecule inhibitor of USP7/HAUSP ubiquitin protease stabilizes and activates p53 in cells. *Mol. Cancer Ther.* 2009; 8:2286–2295. [PubMed: 19671755]
49. Pulvino M, et al. Inhibition of proliferation and survival of diffuse large B-cell lymphoma cells by a small-molecule inhibitor of the ubiquitin-conjugating enzyme Ubc13-Uev1A. *Blood*. 2012; 120:1668–1677. [PubMed: 22791293]
50. Strickson S, et al. The anti-inflammatory drug BAY 11-7082 suppresses the MyD88-dependent signalling network by targeting the ubiquitin system. *The Biochemical journal*. 2013; 451:427–437. [PubMed: 23441730]
51. Daviet L, Colland F. Targeting ubiquitin specific proteases for drug discovery. *Biochimie*. 2008; 90:270–283. [PubMed: 17961905]

52. Reverdy C, et al. Discovery of specific inhibitors of human USP7/HAUSP deubiquitinating enzyme. *Chem. Biol.* 2012; 19:467–477. [PubMed: 22520753]
53. Anderson NL, et al. Precision of heavy-light peptide ratios measured by maldi-tof mass spectrometry. *J. Proteome Res.* 2012; 11:1868–1878. [PubMed: 22257466]
54. Komander D. Mechanism, specificity and structure of the deubiquitinases. *Subcell. Biochem.* 2010; 54:69–87. [PubMed: 21222274]
55. Messick TE, et al. Structural basis for ubiquitin recognition by the Otu1 ovarian tumor domain protein. *The Journal of biological chemistry.* 2008; 283:11038–11049. [PubMed: 18270205]
56. Huang OW, et al. Phosphorylation-dependent activity of the deubiquitinase DUBA. *Nat. Struct. Mol. Biol.* 2012; 19:171–175. [PubMed: 22245969]
57. Virdee S, Ye Y, Nguyen DP, Komander D, Chin JW. Engineered diubiquitin synthesis reveals Lys29-isopeptide specificity of an OTU deubiquitinase. *Nat. Chem. Biol.* 2010; 6:750–757. [PubMed: 20802491]
58. El Oualid F, et al. Chemical synthesis of ubiquitin, ubiquitin-based probes, and diubiquitin. *Angew. Chem. Int. Ed. Engl.* 2010; 49:10149–10153. [PubMed: 21117055]
59. Kramer HB, Nicholson B, Kessler BM, Altun M. Detection of ubiquitin-proteasome enzymatic activities in cells: application of activity-based probes to inhibitor development. *Biochim. Biophys. Acta.* 2012; 1823:2029–2037. [PubMed: 22613766]
60. Pierce JW, et al. Novel inhibitors of cytokine-induced IkappaBalpha phosphorylation and endothelial cell adhesion molecule expression show anti-inflammatory effects in vivo. *The Journal of biological chemistry.* 1997; 272:21096–21103. [PubMed: 9261113]
61. Krishnan N, Bencze G, Cohen P, Tonks NK. The anti-inflammatory compound BAY-11-7082 is a potent inhibitor of protein tyrosine phosphatases. *The FEBS journal.* 2013; 280:2830–2841. [PubMed: 23578302]
62. Lee DW, et al. The Dac-tag, an affinity tag based on penicillin-binding protein 5. *Anal. Biochem.* 2012; 428:64–72. [PubMed: 22705378]
63. Sherrod SD, et al. Label-free quantitation of protein modifications by pseudo selected reaction monitoring with internal reference peptides. *J Proteome Res.* 2012; 11:3467–3479. [PubMed: 22559222]


a. MALDI-TOF DUB Assay workflow

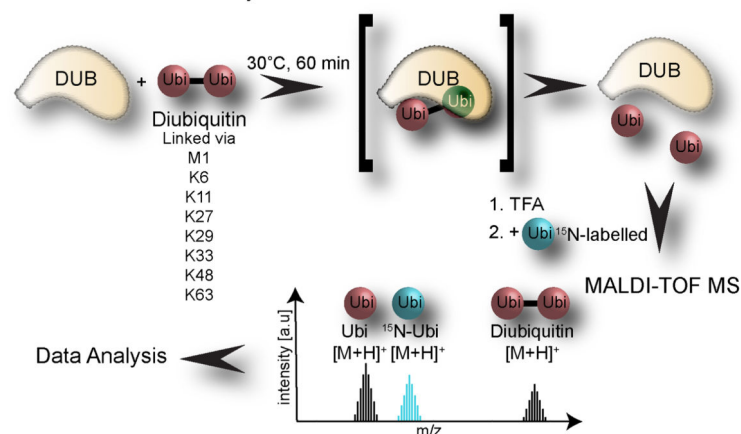

b. Sensitivity

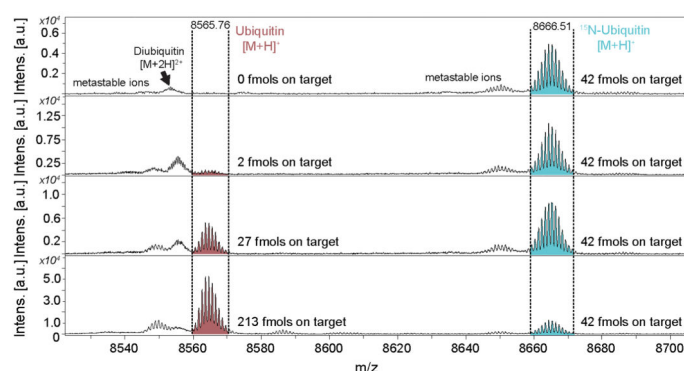

c. Linearity and reproducibility

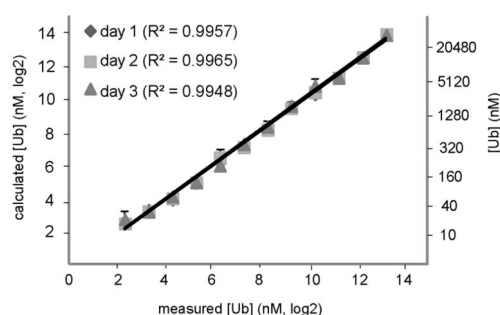**Figure 1. The MALDI-TOF DUB assay**

(a) Workflow of the MALDI-TOF DUB assay. Each of the 42 DUBs was incubated with all eight diubiquitin isomers individually (M1, K6, K11, K27, K29, K33, K48 and K63) for 60 min at 30°C. The reaction was stopped with 2% TFA and mixed 1:1 with 0.5 μ M ¹⁵N-ubiquitin which serves as an internal standard. Subsequently, the analyte is mixed with 2,5 DHAP matrix and spotted onto a 1,536 AnchorChip MALDI target (Bruker Daltonics). Data analysis is performed using FlexAnalysis (Bruker Daltonics). (b) The MALDI-TOF DUB assay shows high sensitivity. Zoomed area (8520–8720 m/z) of MALDI-TOF MS spectra for

ubiquitin and ^{15}N -ubiquitin, in the presence of K11-linked diubiquitin are depicted. The limit of detection was determined as 2 fmol of Ubiquitin on the target (in the presence of 42 fmol of ^{15}N -ubiquitin and 146 fmol of K11-linked diubiquitin). Presence of the doubly-charged diubiquitin (Diubiquitin $[\text{M}+2\text{H}]^{2+}$) does not compromise identification of the singly-charged ubiquitin (**see also** Supplementary Figure 2). (c) Linearity and reproducibility of the MALDI-TOF DUB assay. Scatter plot of different concentrations of ubiquitin (10-10000 nM) shows high linearity over about three orders of magnitude. Interday reproducibility was very high (Supplementary Table 1). Error bars represent standard deviation (SD) of measurements.

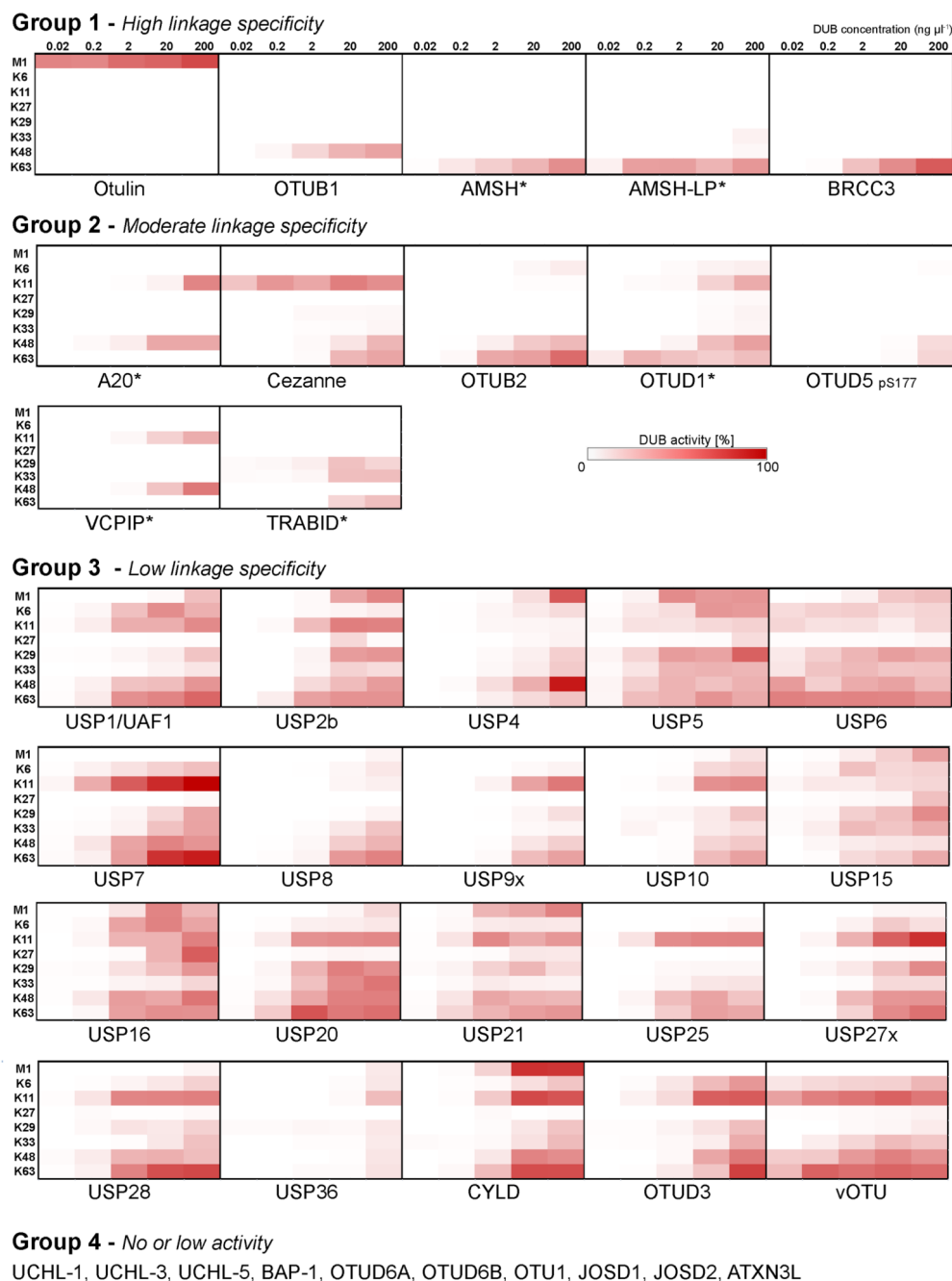


Figure 2. Characterizing the linkage specificity of DUBs

Increasing concentrations (0.02-200 ng μl^{-1}) of DUBs were incubated in triplicate with 1.46 μM of diubiquitin of each linkage type (M1, K6, K11, K29, K33, K48, K63 from Boston Biochem, K27 in-house produced) for 60 min at 30°C and analyzed by the MALDI-TOF DUB assay. The amount of monoubiquitin formed by this reaction was determined by MALDI-TOF MS and used to establish the DUB activity for individual diubiquitin isomers which is shown in a gradient of white (0%) to dark red (100%). The data shows that DUBs can be grouped into enzymes cleaving specifically one linkage type (**group 1**), few linkage

types (**group 2**), unspecific (**group 3**) or inactive enzymes (**group 4**). For DUB characterization, see Supplementary Figures 3 & 8.



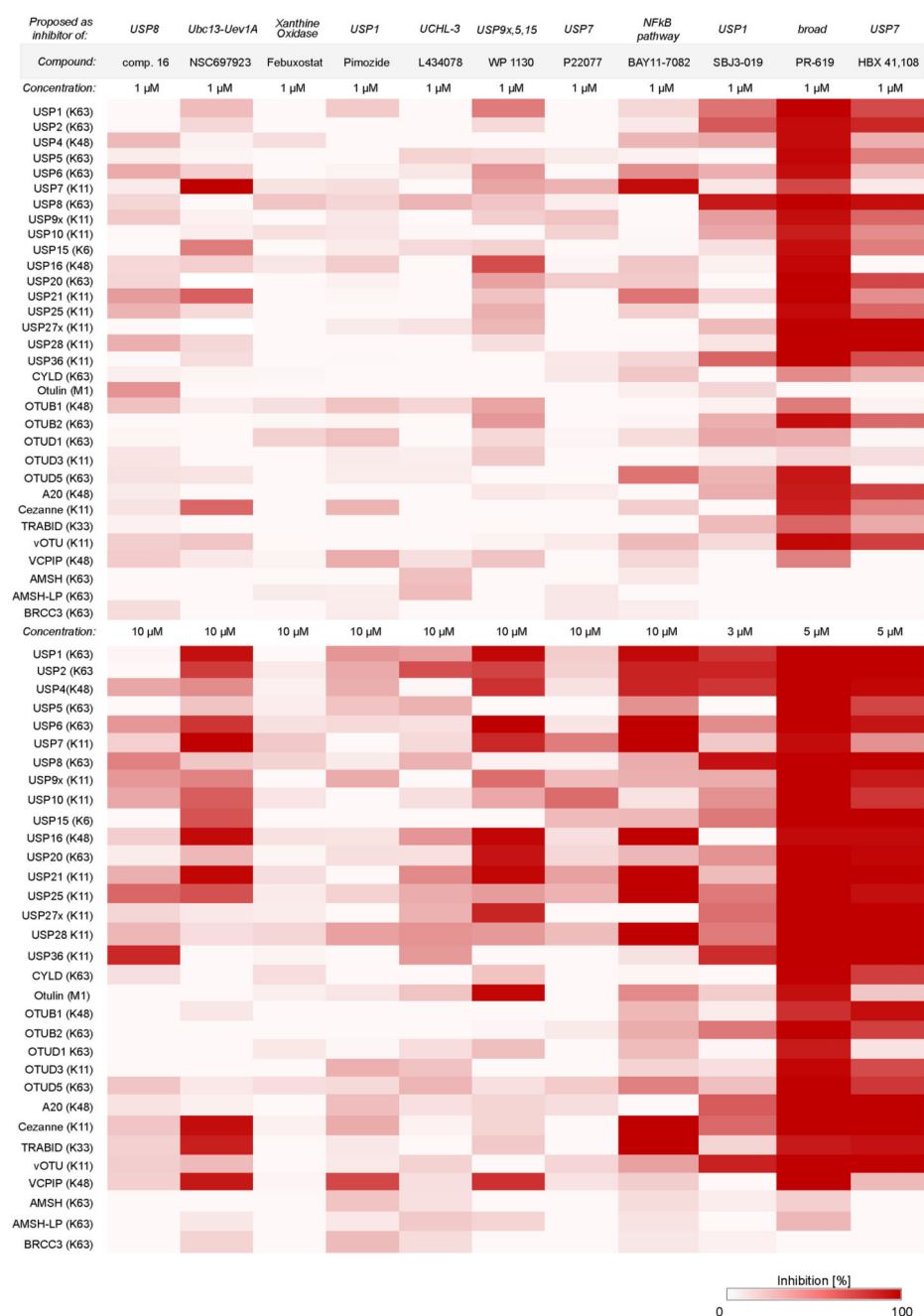


Figure 3. Inhibition profiles of 11 DUB inhibitors and inhibitor candidates

Eleven different DUB inhibitors and inhibitor candidates were pre-incubated for 35 min at two different concentrations in duplicate (i.e. two different experiments) with a panel of 32 DUBs and subsequently the specific substrate was added and incubated for 60 min (30°C). Inhibition rates are color-coded with strongest inhibition in dark red, the diubiquitin topoisomers used for each DUB are in brackets. BAY 11-7082, NSC 697923 and SJB3-019A show some selectivity at 1 μM against USP7 and USP8, respectively, while PR-619 and HBX 41,108 inhibit strongly a wide range of DUBs even at low concentration.

Other proposed inhibitors such as compound 16, L434078, WP1130 and P22077 show low activity and selectivity in this panel.

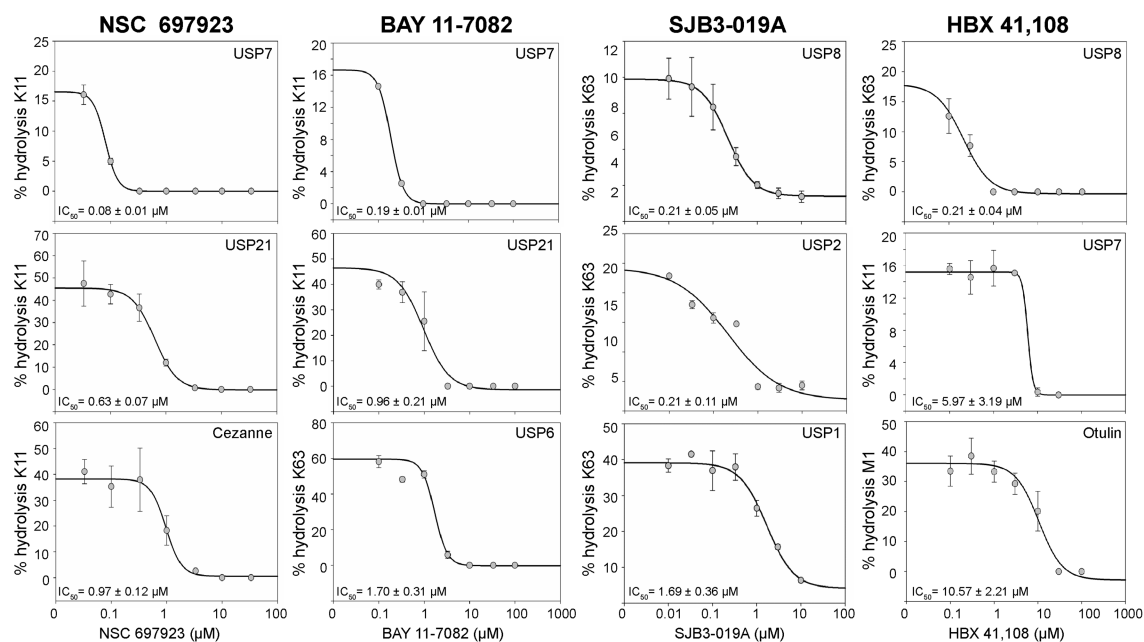


Figure 4. IC₅₀ analyses of four inhibitors for selected DUBs

A subset of four inhibitors was chosen to characterize in more detail by determining their IC₅₀ for three DUBs. BAY 11-7082, NSC 697923 and SJB3-019A were chosen as they have some selectivity for one DUB, HBX 41,108 was chosen as it has been proposed as a USP7 inhibitor which is an attractive drug target⁵¹. Small inhibitor compounds were pre-incubated for 35 min at different concentrations (0–30 μM or 0–100 μM) in triplicates (i.e. three different experiments) and subsequently the specific substrate was added and incubated for 60 min (30°C). Diubiquitin topoisomers used for each DUB are named on the y-axis. Data shows that NSC 697923 and BAY 11-7082 inhibit strongly USP7 with IC₅₀ < 0.2 μM, while HBX 41,108 inhibits it at ~6 μM. SJB3-019A inhibits USP8 and USP2 about ten-fold better than USP1. See also Supplementary Table 4 for p-values. Error bars represent standard deviation (SD) of measurements. For statistical analysis, four parameter logistic curve (Best - Fit Solution, Nonlinear Regression-Dynamic Fitting) and normality tests (Kolmogorov-Smirnov) are used (SigmaPlot, v. 12.5).

Table 1
DUB enzymes analyzed in this study

	DUB family	DUB	UniProt Accession Number	Tag	Domain	Host	MALDI-TOF DUB	Rhodamine-110-Gly	
							DUB concentration for hydrolysis rate≥20%	Specific activity (counts/min/ng)	Activity ¹
1	ubiquitin- specific proteases (USPs)	USP1/UAF1	O94782	His	Full length	<i>S. frugiperda</i>	2 ng μl ⁻¹	92	**
2		USP2b	O75604	GST	Full length	<i>E. coli</i>	2 ng μl ⁻¹	133	***
3		USP4	Q13107	His	Full length	<i>E. coli</i>	20 ng μl ⁻¹	367	***
4		USP5	P45974	His	Full length	<i>E. coli</i>	0.2 ng μl ⁻¹	59	**
5		USP6	P35125	GST	[529 ⁻ 1406]	<i>S. frugiperda</i>	0.02 ng μl ⁻¹	638	***
6		USP7	Q93009	His	Full length	<i>E. coli</i>	0.2 ng μl ⁻¹	15	**
7		USP8	P40818	His	Full length	<i>S. frugiperda</i>	20 ng μl ⁻¹	16	**
8		USP9x	Q93008	GST	[1553-1995]	<i>E. coli</i>	20 ng μl ⁻¹	6	-
9		USP10	Q14694	His	Full length	<i>S. frugiperda</i>	20 ng μl ⁻¹	10	*
10		USP15	Q9Y4E8	GST	Full length	<i>E. coli</i>	2 ng μl ⁻¹	107	***
11		USP16	Q9Y5T5	His	Full length	<i>E. coli</i>	2 ng μl ⁻¹	323	***
12		USP20	Q9Y2K6	GST	Full length	<i>S. frugiperda</i>	2 ng μl ⁻¹	405	***
13		USP21	Q9UK80	His	[196-565]	<i>E. coli</i>	2 ng μl ⁻¹	416	***
14		USP25	Q9UHP3	GST	Full length	<i>E. coli</i>	2 ng μl ⁻¹	17	**
15		USP27x	A6NNY8	DAC	Full length	<i>S. frugiperda</i>	2 ng μl ⁻¹	8	-
16		USP28	Q96RU2	GST	Full length	<i>E. coli</i>	2 ng μl ⁻¹	10	*
17		USP36	Q9P275	GST	[81-461]	<i>E. coli</i>	200 ng μl ⁻¹	6	-
18		CYLD	Q9NQC7	His	Full length	<i>S. frugiperda</i>	2 ng μl ⁻¹	4	-
19	UCHs	UCHL1	P09936	His	Full length	<i>E. coli</i>	-	218	***
20		UCHL3	P15374	GST	Full length	<i>E. coli</i>	-	623	***
21		UCHL5	Q9Y5K5	GST	Full length	<i>E. coli</i>		256	***
22		BAP1	Q92560	GST	Full length	<i>E. coli</i>		46	**
23	ovarian tumor proteases (OTUs)	OTULIN	Q96BN8	GST	Full length	<i>E. coli</i>	0.02 ng μl ⁻¹	4	-
24		OTUB1	Q96FW1	GST	Full length	<i>E. coli</i>	20 ng μl ⁻¹	3	-
25		OTUB2	Q96DC9	GST	Full length	<i>E. coli</i>	2 ng μl ⁻¹	5	-
26		OTUD1	Q5VV17	His	[270-481]	<i>E. coli</i>	0.2 ng μl ⁻¹	73	**
27		OTUD3	Q5T2D3	His	Full length	<i>E. coli</i>	20 ng μl ⁻¹	255	***
28		OTUD5 pS177	Q96G74	GST	Full length	<i>E. coli</i>	200 ng μl ⁻¹	18	**
29		OTUD6A	Q7L8S5	His	Full length	<i>E. coli</i>	-	539	***
30		OTUD6B	Q8N6M0	GST	Full length	<i>E. coli</i>	-	9	-
31		OTU1	Q5VVQ6	GST	Full length	<i>E. coli</i>	-	8	-
32		A20	P21580	GST	[1-366]	<i>E. coli</i>	20 ng μl ⁻¹	9	*
33		Cezanne	Q6GQQ9	GST	Full length	<i>E. coli</i>	0.02 ng μl ⁻¹	466	***

DUB family	DUB	UniProt Accession Number	Tag	Domain	Host	MALDI-TOF DUB	Rhodamine-110-Gly	Activity ¹
						DUB concentration for hydrolysis rate ≥ 20%	Specific activity (counts/min/ng)	
34	TRABID	Q9UGI0	His	[245-697]	<i>E. coli</i>	20 ng μl^{-1}	16	**
35	vOTU	Q6TQR6	His	[1-183]	<i>E. coli</i>	0.02 ng μl^{-1}	670	***
36	VCPIP1	Q96JH7	His	Full length	<i>E. coli</i>	2 ng μl^{-1}	13	*
36	VCPIP1	Q96JH7	GST	[25-561]	<i>E. coli</i>	20 ng μl^{-1}	12	*
37	JOSD1	Q15040	His	Full length	<i>E. coli</i>	-	20	**
38	MJD+ JOSD2	Q8TAC2	His	Full length	<i>E. coli</i>	-	175	***
39	ATXN3L	Q9H3M9	His	Full length	<i>E. coli</i>	-	8	-
40	AMSH	O95630	GST	[256-424]	<i>E. coli</i>	20 ng μl^{-1}	5	-
41	JAMM/MPN+ AMSH-LP	Q96FJ0	GST	[265-436]	<i>E. coli</i>	0.2 ng μl^{-1}	6	-
42	BRCC3	P46736	His	Full length	<i>E. coli</i>	2 ng μl^{-1}	8	-

¹(-) very low (0-8.9 counts $\text{min}^{-1} \text{ng}^{-1}$) (*) low (9-14 counts $\text{min}^{-1} \text{ng}^{-1}$) (**) moderate (15-99 counts $\text{min}^{-1} \text{ng}^{-1}$) (***) high (100-700 counts $\text{min}^{-1} \text{ng}^{-1}$) activity by Rhodamine-110-Glycine assay.